


InertSearch for LC Data No. LA272-0000

Pesticides

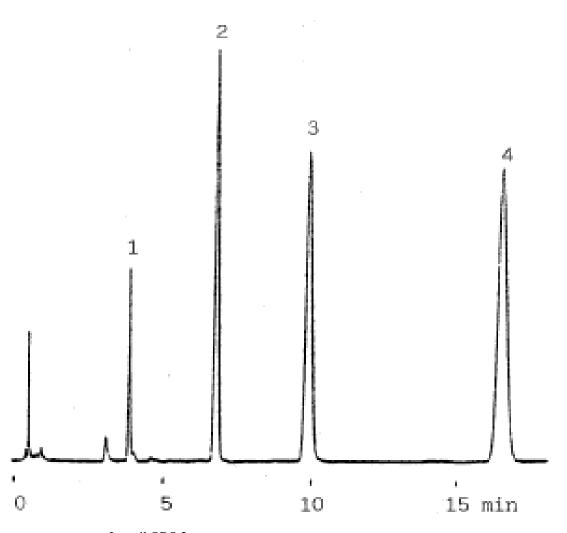

Column

: Inertsil ODS-2

5µm

150 × 4.6mmI.D.

Column Cat. No.: 5020-01101

Eluent

: CH3CN/H2O=65/35

Flow rate

: 1.0 mL/min.

Col. Temp.

:40°C

Detection

:UV 210 nm

Injection Vol.

Sample

Solutes

: 1. Fenobucarb(BPMC)

2. Bensulide

3. Isofenphos

4. Pendimethalin