


# InertSearch™ for LC

Inertsil® Applications

## Pesticides

Data No. LA274-0000

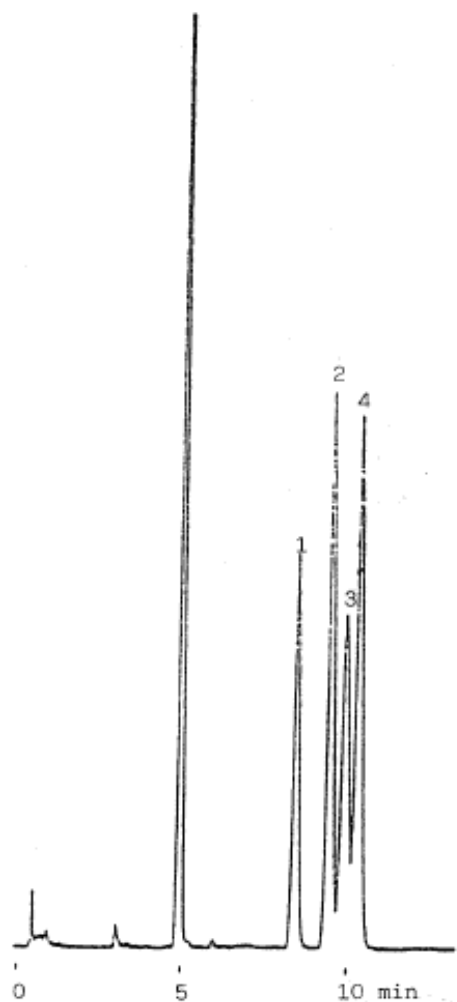

Column : Inertsil ODS-2  
5µm  
150 × 4.6mm I.D.

Column Cat. No. : 5020-01101

Eluent : CH<sub>3</sub>CN/H<sub>2</sub>O=65/35

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 210 nm

Injection Vol. :

Sample :

Solutes : 1. Flutolanil  
2. Isoxathion  
3. Diazinon  
4. Butamifos