


InertSearch™ for LC

Inertsil® Applications

Pesticides

Data No. LA276-0000

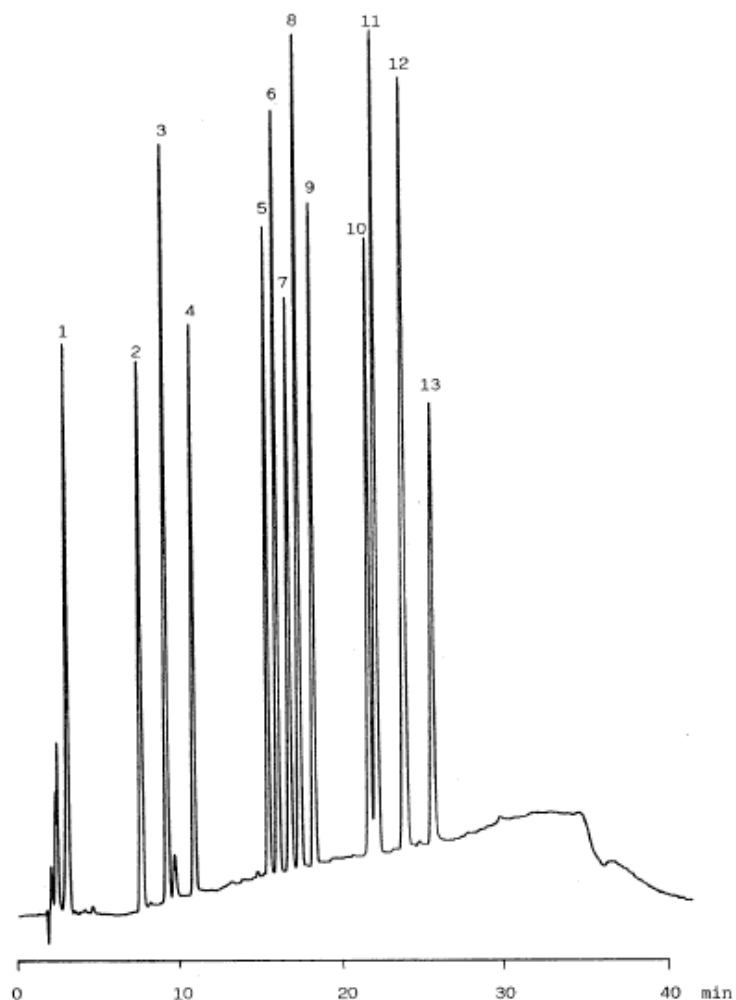

Column	: Inertsil C4 5µm 150 × 4.6mm I.D.	Solutes	: 1. Asulam 20µg/mL 2. Simazine 20µg/mL 3. DDVP 200µg/mL 4. Thiuram 20µg/mL 5. Fenobucarb 300µg/mL 6. Captan 200µg/mL 7. Propyzamide 20µg/mL 8. Napropamide 10µg/mL 9. Iprodione 20µg/mL 10. Tolclophos-methyl 40µg/mL 11. Bensulide 70µg/mL 12. Butamifos 80µg/mL 13. Chlorpyrifos 50µg/mL
Column Cat. No.	: 5020-01401		
Eluent	: A: CH ₃ CN B: 0.1% H ₃ PO ₄ A/B 25/75 → (25min) → 75/25		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 210 nm		
Injection Vol.	:		
Sample	:		