


InertSearch™ for LC

Inertsil® Applications

Pesticides

Data No. LA278-0000

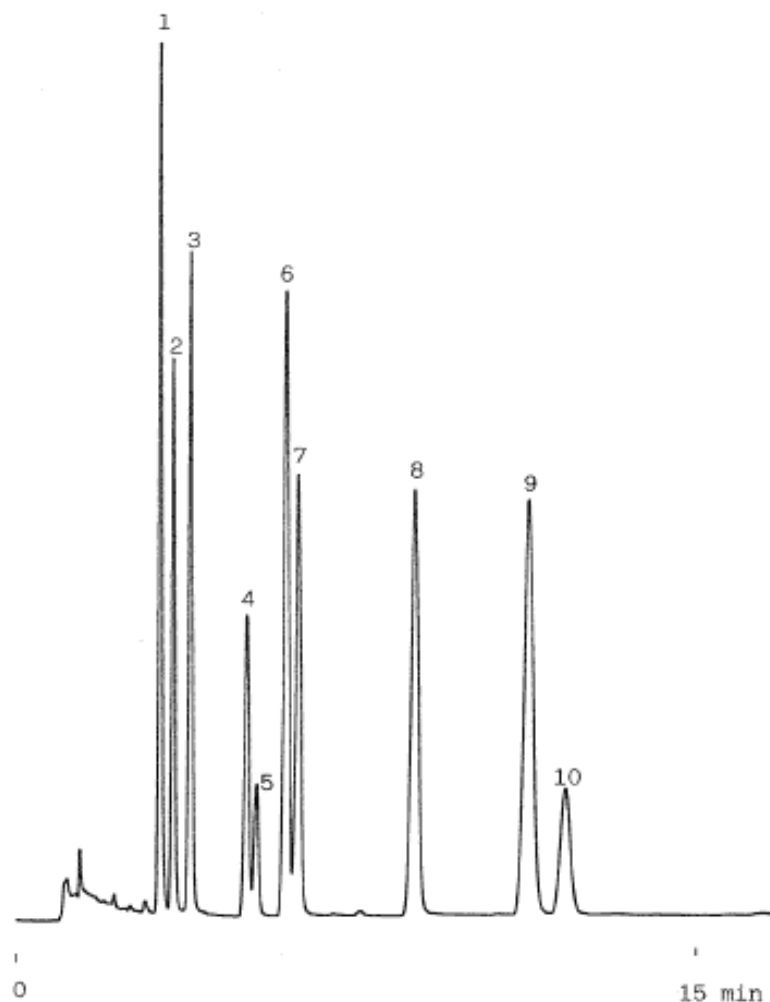

Column : Inertsil ODS-2
5 μ m
150 \times 4.6mm I.D.

Column Cat. No. : 5020-01101

Eluent : CH₃CN/0.1% H₃PO₄=60/40

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 230 nm

Injection Vol. :

Sample :

Solutes : 1. MCPA
2. Triclopyr
3. MCPP
4. Pyridaphenthion
5. Methyl Dymron
6. Mepronil
7. Chloroneb
8. Echlomezol
9. Pencycuron
10. Terbutcarb