


# Pesticides

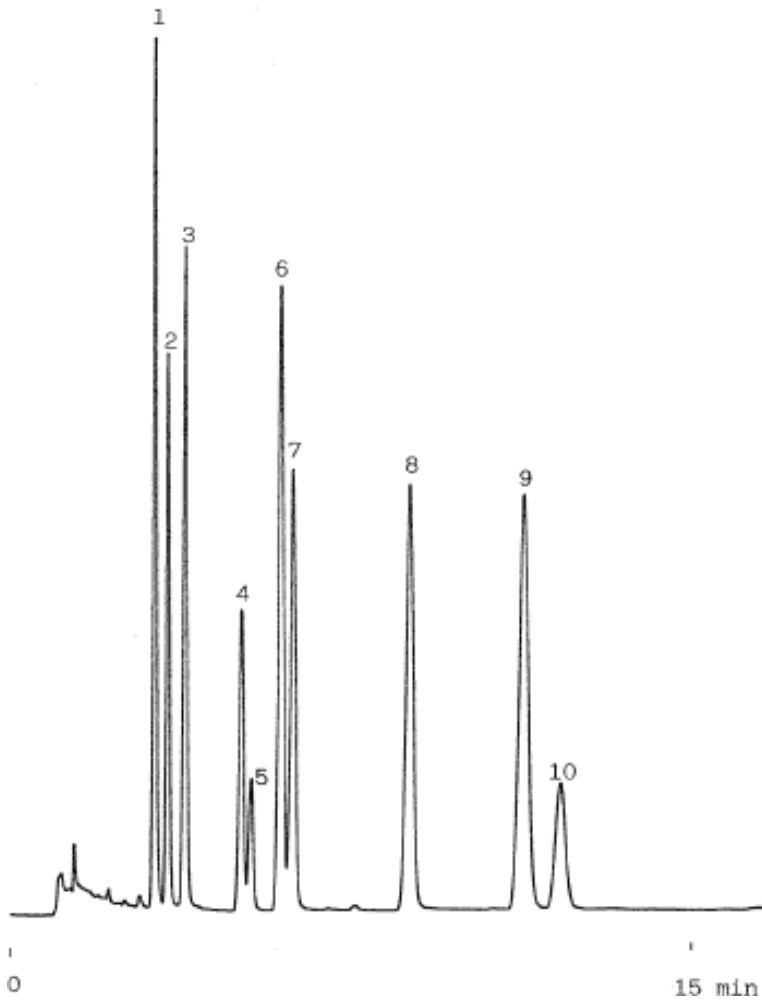

Column	: Inertsil ODS-2 5µm 150 × 4.6mmI.D.	Solutes	: 1. MCPA 2. Triclopyr 3. MCPP 4. Pyridaphenthion 5. Methyl Dymron 6. Mepronil 7. Chloroneb 8. Echlomezol 9. Pencycuron 10. Terbutcarb
Column Cat. No.	: 5020-01101		
Eluent	: CH <sub>3</sub> CN/0.1% H <sub>3</sub> PO <sub>4</sub> =60/40		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 230 nm		
Injection Vol.	:		
Sample	:		