


Pesticides

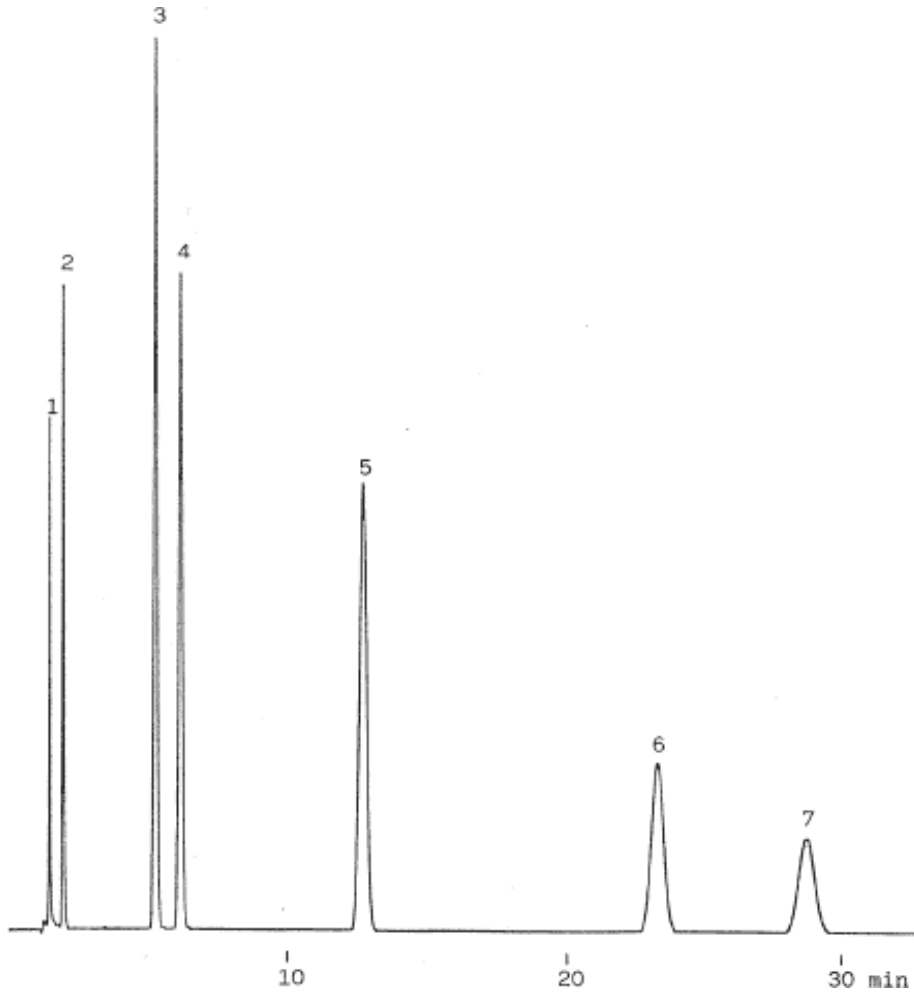

Column : Inertsil ODS-3
 5µm
 150 × 4.6mmI.D.

Column Cat. No. : 5020-13001

Eluent : CH₃CN/50mM KH₂PO₄(pH2.0)=50/50

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 230 nm

Injection Vol. :

Sample :

Solutes : 1. Asulam
 2. Oxine-copper
 3. Thiuram
 4. MCP
 5. Iprodione
 6. Bensulide
 7. Pencycuron