


InertSearch™ for LC

Inertsil® Applications

Pesticides

Data No. LA289-0000

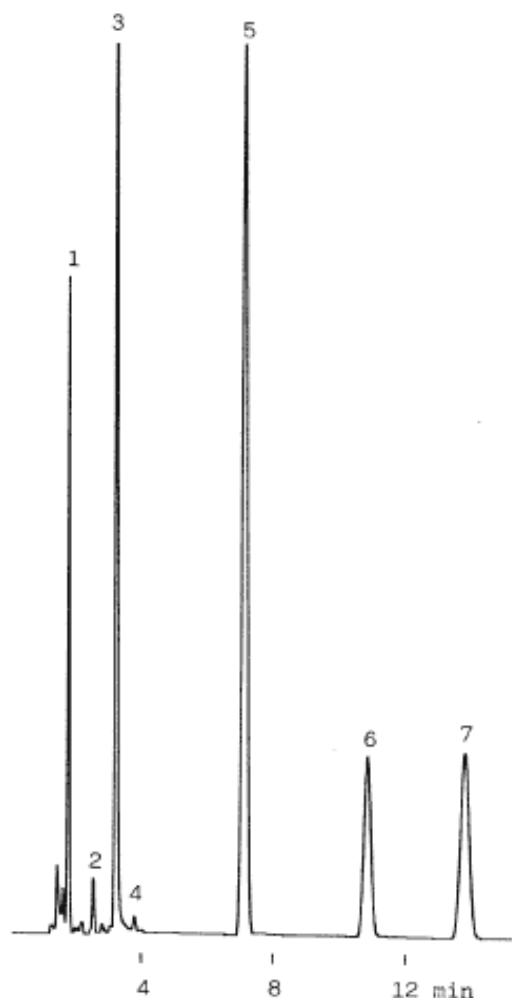

Column : Inertsil ODS-3
5 μ m
150 \times 4.6mm I.D.

Column Cat. No. : 5020-13001

Eluent : CH₃CN/50mM KH₂PO₄(pH4.0 by H₃PO₄)=60/40

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 230 nm

Injection Vol. :

Sample :

Solutes : 1. Asulam
2. MCPP
3. Oxine-copper
4. Thiuram
5. Iprodione
6. Bensulide
7. Penycuron