


# Pesticides

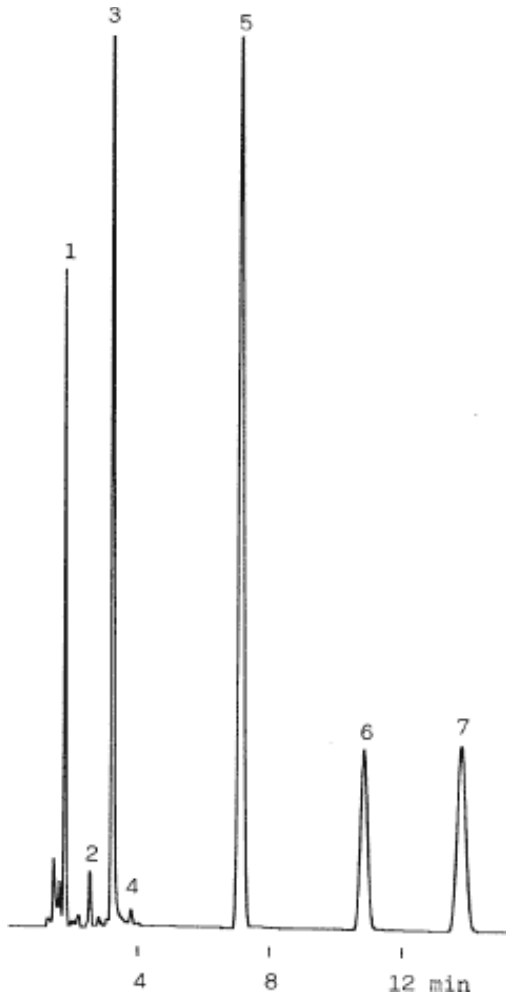

Column : Inertsil ODS-3  
 5µm  
 150 × 4.6mmI.D.

Column Cat. No. : 5020-13001

Eluent : CH<sub>3</sub>CN/50mM KH<sub>2</sub>PO<sub>4</sub>(pH4.0 by H<sub>3</sub>PO<sub>4</sub>)=60/40

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 230 nm

Injection Vol. :

Sample :

Solutes : 1. Asulam  
 2. MCPP  
 3. Oxine-copper  
 4. Thiuram  
 5. Iprodione  
 6. Bensulide  
 7. Pencycuron