


Analysis of Pesticides

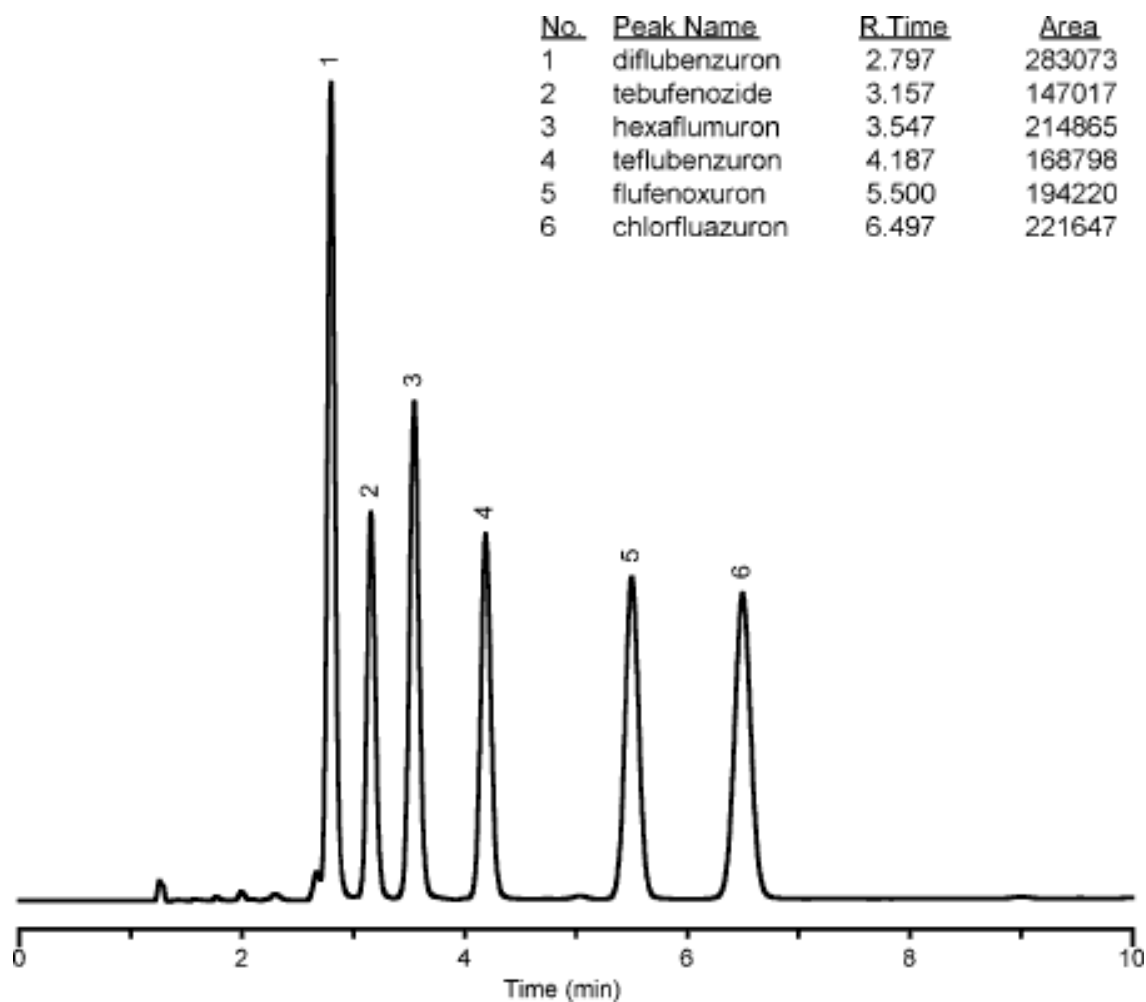

Column	: Inertsil ODS-3V 5µm 150 × 4.6mmI.D.	Solutes	: 1. Diflubenzuron (10 g/mL) 2. Tebufenozide (10 g/mL) 3. Hexaflumuron (10 g/mL) 4. Teflubenzuron (10 g/mL) 5. Flufenoxuron (10 g/mL) 6. Chlorfluazuron (10 g/mL)
Column Cat. No.	: 5020-01801		
Eluent	: Acetonitrile / Water = 80 / 20		
Flow rate	: 1mL/min		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	: 10µL		
Sample	: Pesticides		