


InertSearch™ for LC

Inertsil® Applications

Analysis of Pesticides

Data No. LA295-0000

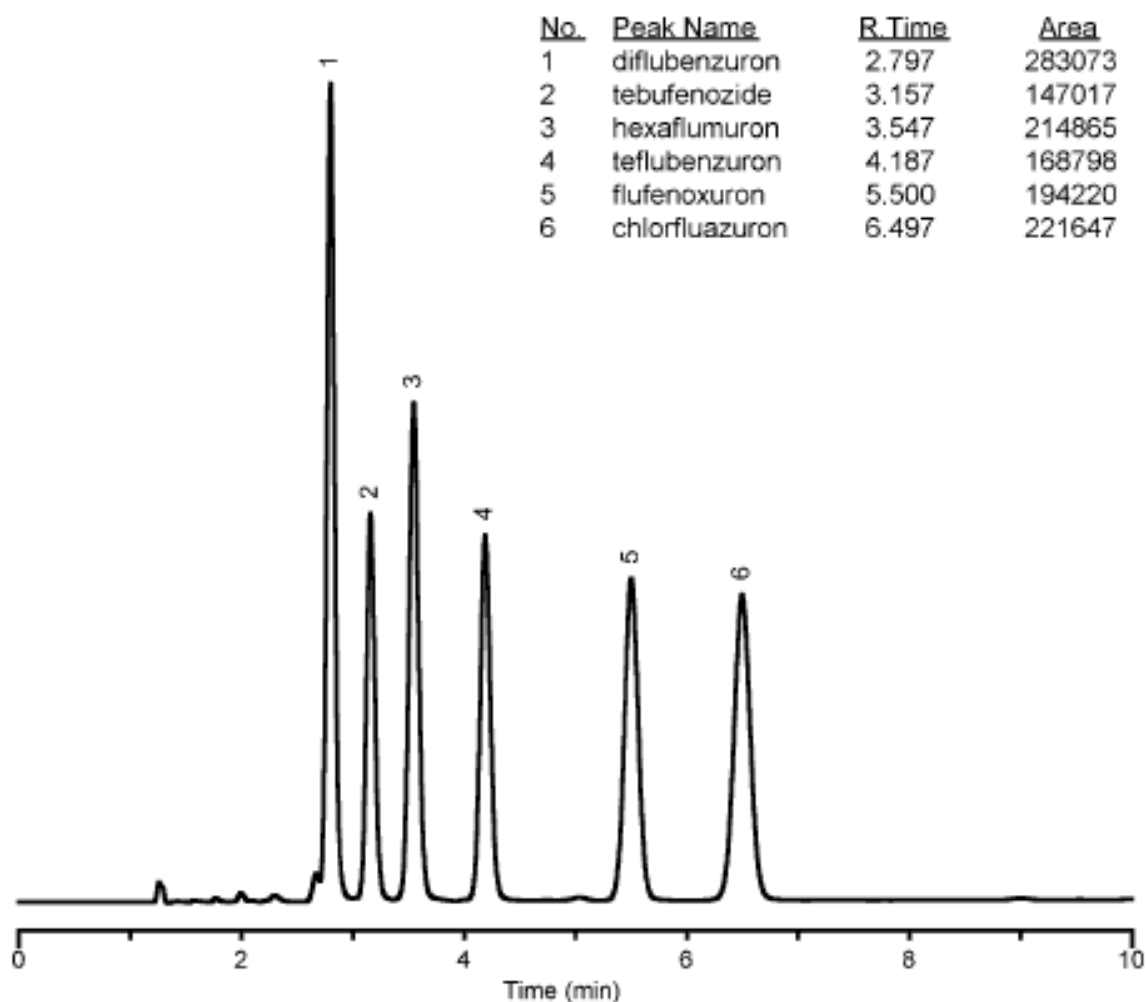

Column : Inertsil ODS-3V
5µm
150 × 4.6mm I.D.

Column Cat. No. : 5020-01801

Eluent : Acetonitrile / Water = 80 / 20

Flow rate : 1mL/min

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. : 10µL

Sample : Pesticides

Solutes : 1. Diflubenzuron (10 g/mL)
2. Tebufenozide (10 g/mL)
3. Hexaflumuron (10 g/mL)
4. Teflubenzuron (10 g/mL)
5. Flufenoxuron (10 g/mL)
6. Chlorfluazuron (10 g/mL)