


## Analysis of Carbmates

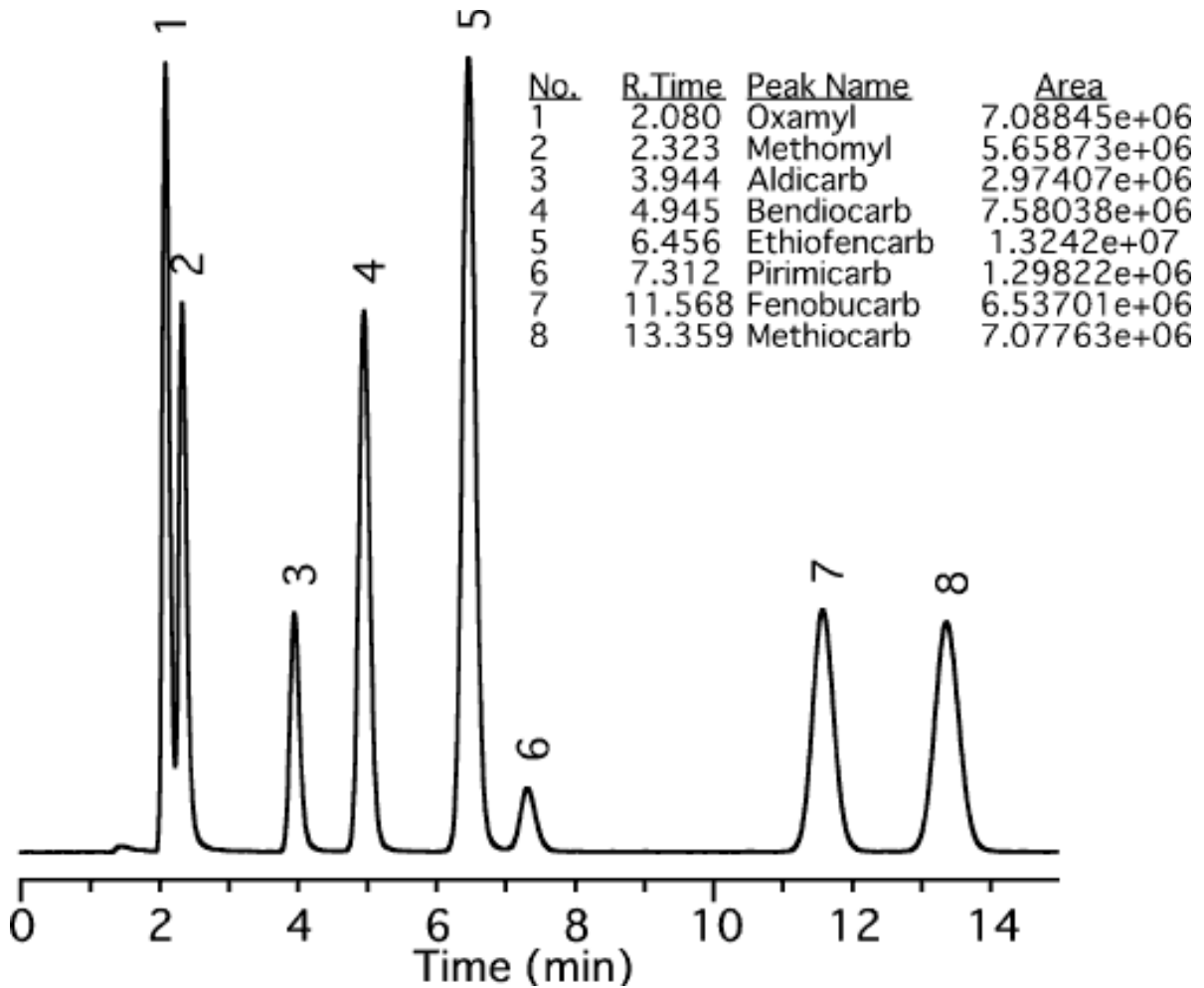

Column : Inertsil ODS-3  
5µm  
150 × 4.6mmI.D.

Column Cat. No. :

Eluent : 53% MeOH

Flow rate : 1mL/min

Col. Temp. : 40°C

Detection : FL Ex 345nm Em 450nm

POST REACTION NaOH 100? to OPA 50?

Injection Vol. : 25µL

Sample : Pesticides

Solutes

- : 1 . Oxamyl 2.83mg
  - 2 . Methomyl 2.5mg
  - 3 . Aldicarb 1.65mg
  - 4 . Bendiocarb 2.5mg
  - 5 . Ethiofencarb 5.10mg
  - 6 . Pirimicarb 7.8mg
  - 7 . Fenobucarb 2.0mg
  - 8 . Methiocarb 2.4mg
- in 50mL 52% MeOH