


InertSearch™ for LC

Inertsil® Applications

Analysis of Oxine-Copper

Data No. LA298-0000

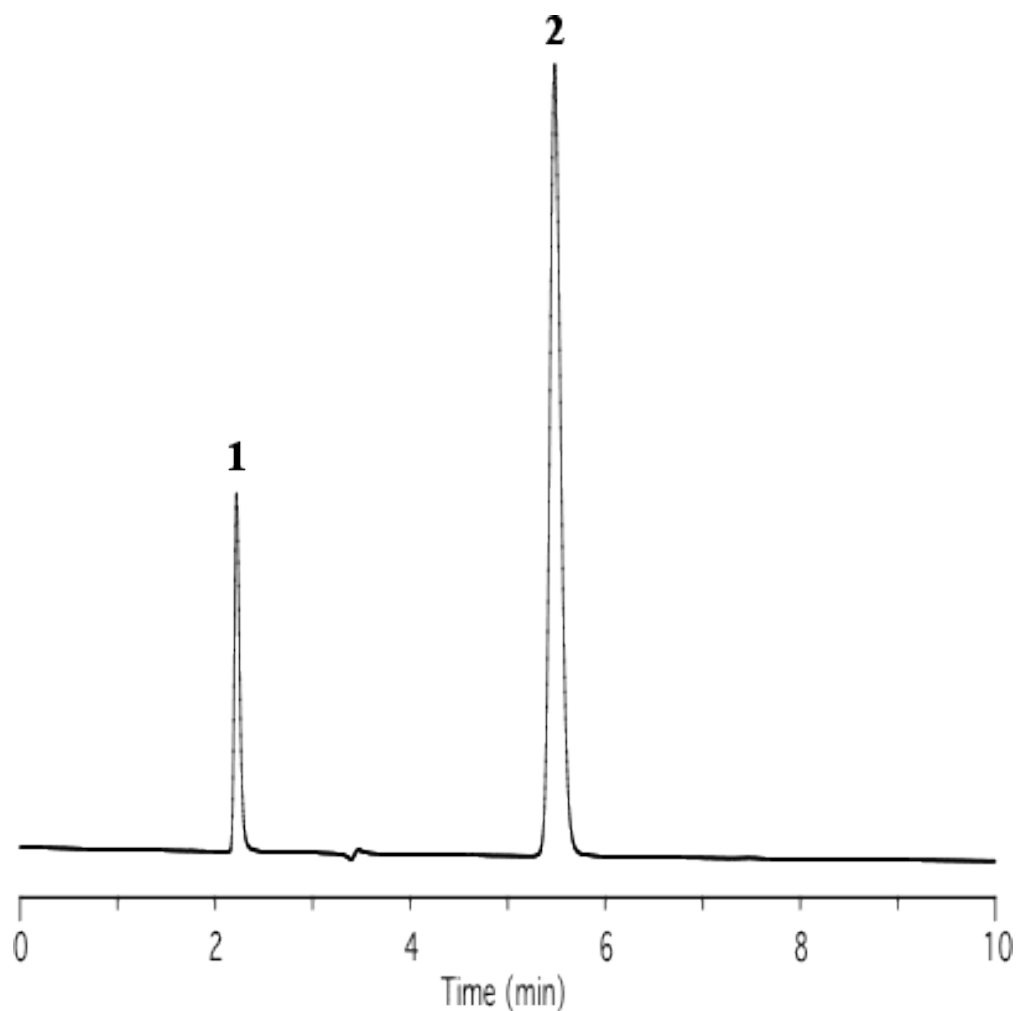

Column : Inertsil CN-3
5 µm
250 x 4.6 mm I.D.

Solutes : 1. Oxine-copper
2. Caffeine

Column Cat. No. :

Eluent : A) CH₃CN
B) 20 mM H₃PO₄
A/B = 5/95, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :