


InertSearch™ for LC

Inertsil® Applications

Analysis of Pesticides

Data No. LA301-0000

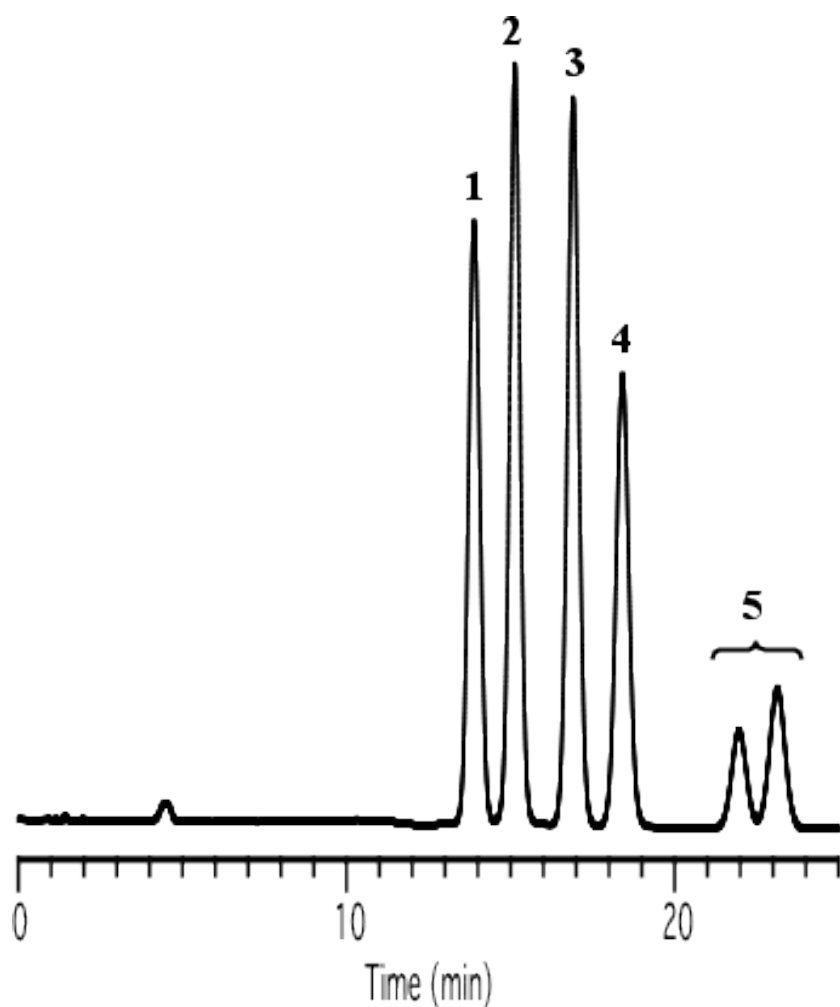

Column	: Inertsil Ph-3 5 µm 150 x 4.6 mm I.D.	Solutes	: 1. Flucythrinate 2. Cyhalothrin 3. Deltamethrin 4. Fluvalinate 5. Tralomethrin
Column Cat. No.	:		
Eluent	: A) CH ₃ CN B) H ₂ O A/B = 60/40, v/v		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 220 nm		
Injection Vol.	:		
Sample	:		