


## Analysis of Pesticides

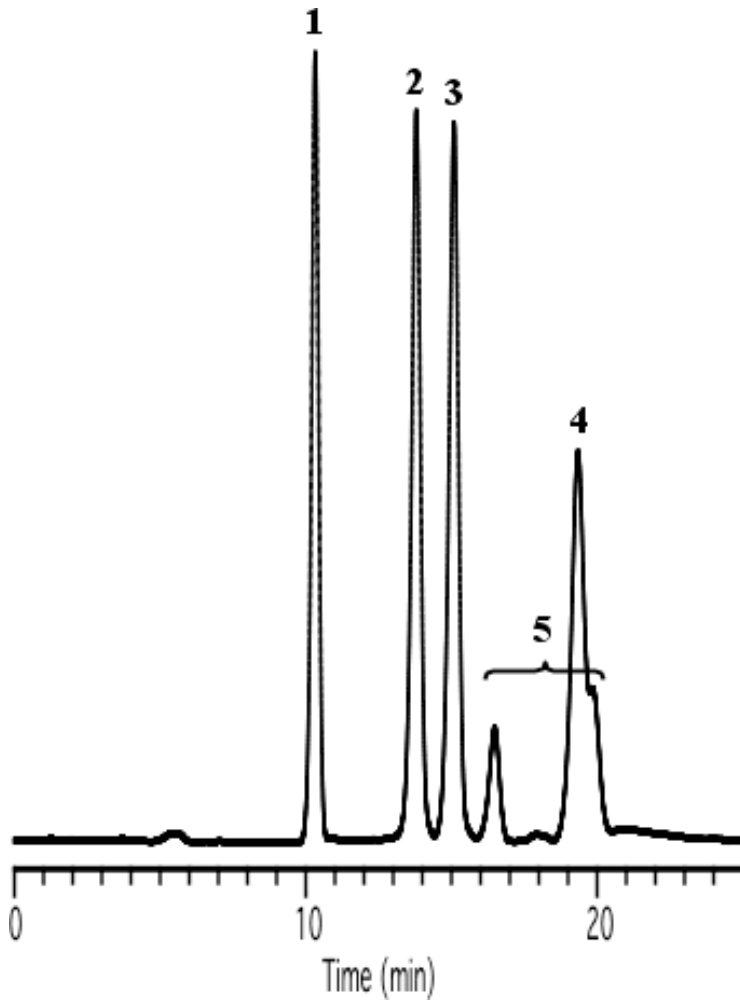

|                   |                                                                    |         |                                                                                              |
|-------------------|--------------------------------------------------------------------|---------|----------------------------------------------------------------------------------------------|
| Column            | : Inertsil ODS-3<br>5 $\mu$ m<br>150 x 4.6 mm I.D.                 | Solutes | : 1. Flucythrinate<br>2. Cyhalothrin<br>3. Deltamethrin<br>4. Fluvalinate<br>5. Tralomethrin |
| Column Cat. No. : |                                                                    |         |                                                                                              |
| Eluent            | : A) CH <sub>3</sub> CN<br>B) H <sub>2</sub> O<br>A/B = 75/25, v/v |         |                                                                                              |
| Flow rate         | : 1.0 mL/min.                                                      |         |                                                                                              |
| Col. Temp.        | : 40°C                                                             |         |                                                                                              |
| Detection         | : UV 220 nm                                                        |         |                                                                                              |
| Injection Vol.    | :                                                                  |         |                                                                                              |
| Sample            | :                                                                  |         |                                                                                              |