


Polyaromatic hydrocarbons

Data No. LA305-0000

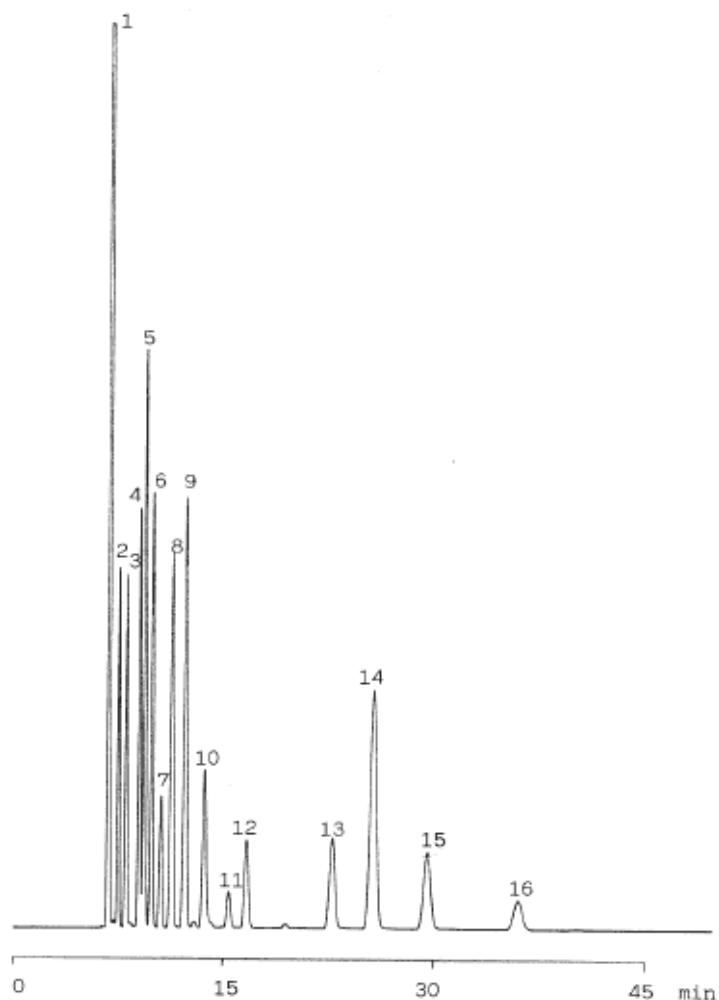

Column : Inertsil ODS-2
5µm
250 × 4.6mm I.D.

Column Cat. No. : 5020-01102

Eluent : CH₃CN/H₂O=75/25

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 220 nm

Injection Vol. :

Sample :

Solutes

- : 1. Naphthalene
- 2. Acenaphthylene
- 3. Biphenyl (Diphenyl)
- 4. Fluorene
- 5. Acenaphthene
- 6. Phenanthrene
- 7. Anthracene
- 8. Benzo(a)anthracene
- 9. Fluoranthene
- 10. Pyrene
- 11. Triphenylene
- 12. Crycene
- 13. Benzo(e)pyrene
- 14. Benzo(a)pyrene
- 15. Dibenzo(ah)anthracene
- 16. Benzo(ghi)perylene