


Chlorobenzene

Data No. LA307-0000

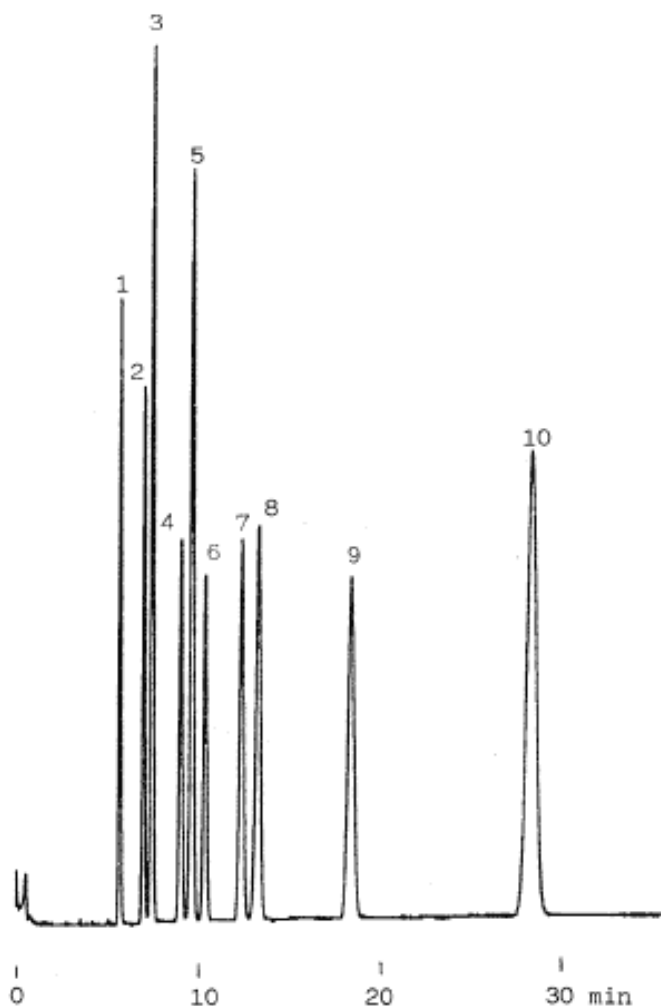

Column : Inertsil Ph
5 μ m
150 \times 4.6mm I.D.

Column Cat. No. : 5020-01301

Eluent : CH₃CN/H₂O=50/50

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes

1. Chlorobenzene
2. 1,2-Dichlorobenzene
3. 1,4-Dichlorobenzene
4. 1,3-Dichlorobenzene
5. 1,2,3-Trichlorobenzene
6. 1,3,5-Trichlorobenzene
7. 1,2,3,4-Tetrachlorobenzene
8. 1,2,4,5-Tetrachlorobenzene
+1,2,3,5-Tetrachlorobenzene
9. Pentachlorobenzene
10. Hexachlorobenzene