


Analysis of DNPH-Aldehydes

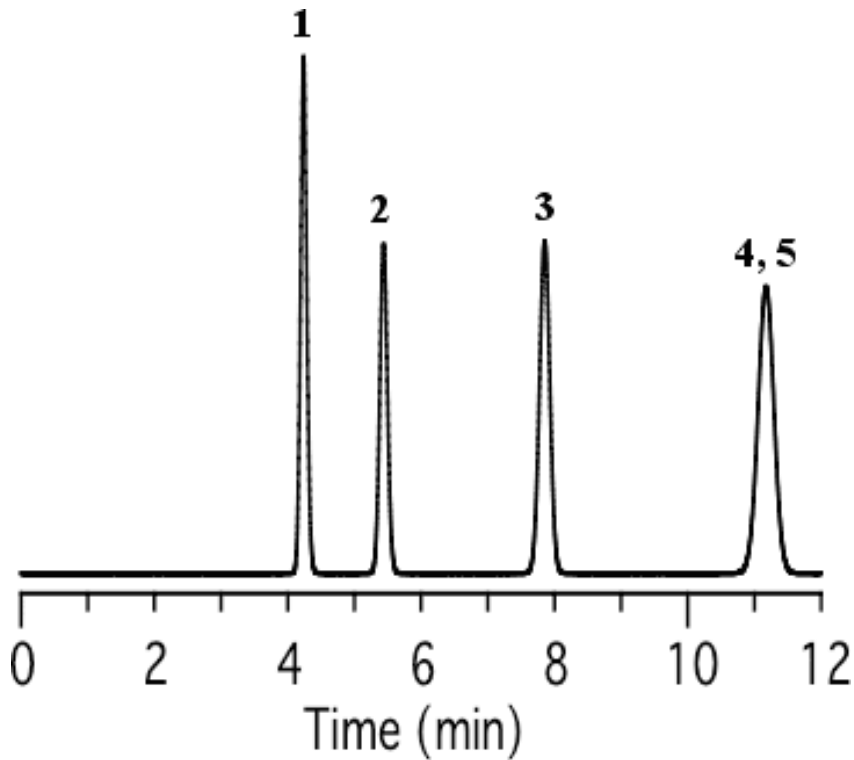

Column	: Inertsil ODS-3V 5 µm 150 x 4.6 mm I.D.	Solutes	: 1. DNPH-Formaldehyde 2. DNPH-Acetaldehyde 3. DNPH-Propionaldehyde 4. DNPH-iso-Butylaldehyde 5. DNPH-n-Butylaldehyde
Column Cat. No. :			
Eluent	: A) CH ₃ CN B) H ₂ O A/B = 60/40, v/v		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 360 nm		
Injection Vol.	:		
Sample	:		