


InertSearch™ for LC

Inertsil® Applications

Analysis of Aromatic Compounds

Data No. LA349-0000

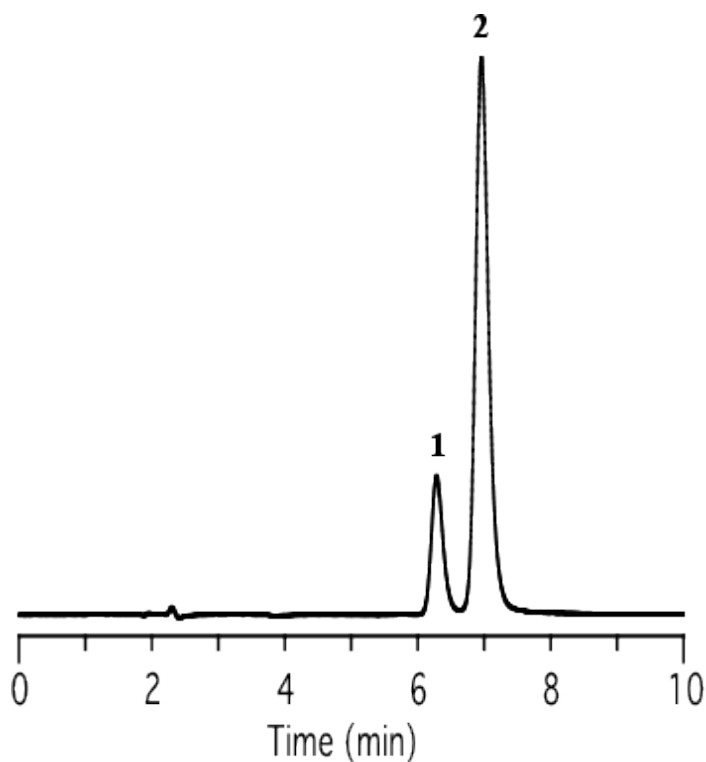

Column : Inertsil Ph-3
5 μ m
150 x 4.6 mm I.D.

Solutes : 1. Fluorene
2. Phenanthrene

Column Cat. No. : 5020-01920

Eluent : A) CH₃OH
B) H₂O
A/B = 70/30, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 30°C

Detection : UV 254 nm

Injection Vol. :

Sample :