


## Cation surfactants

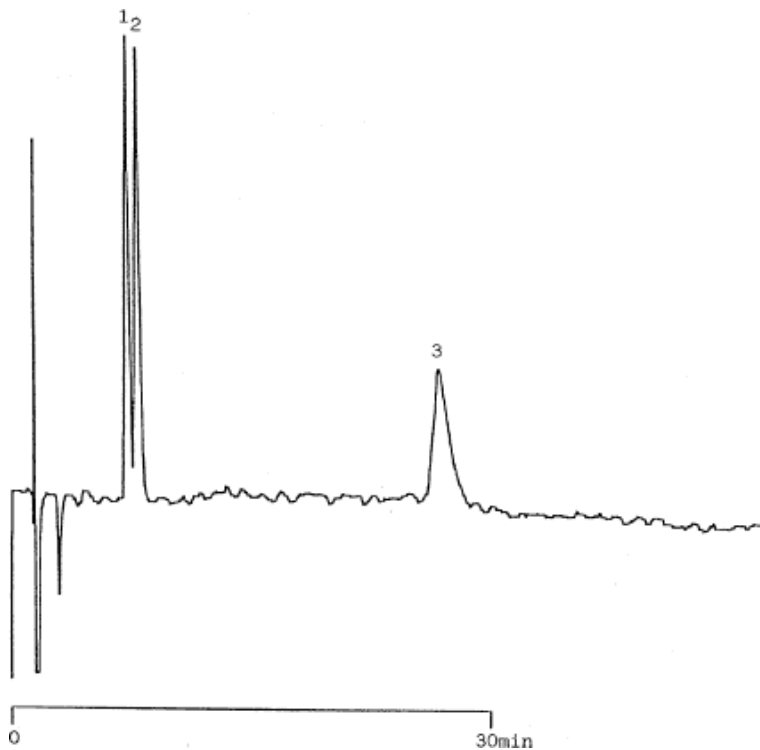

Column	: Inertsil C8 5 $\mu$ m 150 $\times$ 4.6 mmI.D.	Solutes	: 1. Undecyl ammonium chloride 2. Lauryl pyridinium chloride 3. Cetyl pyridinium chloride
Column Cat. No.	: 5020-01201		
Eluent	: CH <sub>3</sub> CN/0.1M NaClO <sub>4</sub> =60/40		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: RI		
Injection Vol.	:		
Sample	:		