


InertSearch™ for LC

Inertsil® Applications

Cation surfactants

Data No. LA352-0000

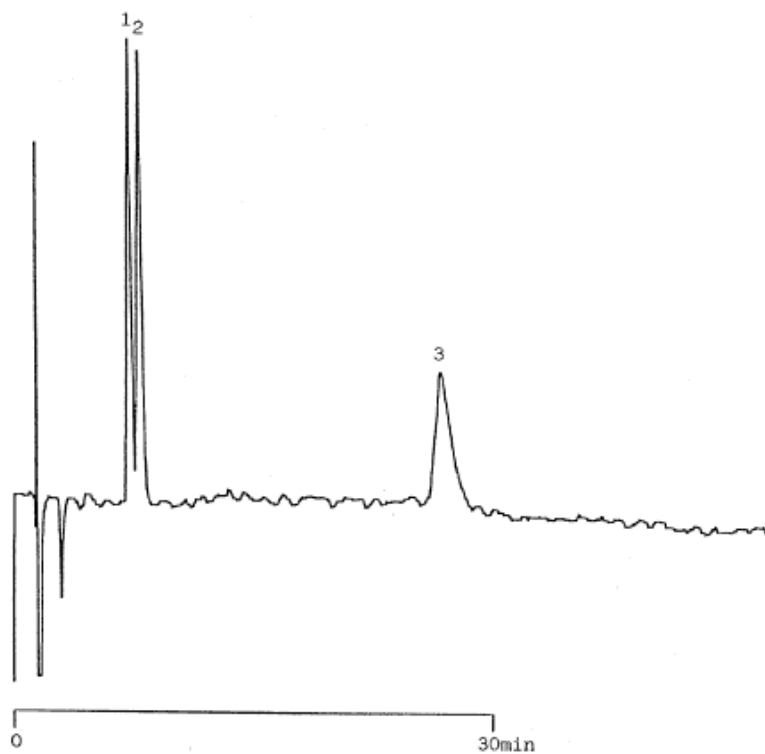

Column	: Inertsil C8 5µm 150 × 4.6 mmI.D.	Solutes	: 1. Undecyl ammonium chloride 2. Lauryl pyridinium chloride 3. Cetyl pyridinium chloride
Column Cat. No.	: 5020-01201		
Eluent	: CH ₃ CN/0.1M NaClO ₄ =60/40		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: RI		
Injection Vol.	:		
Sample	:		