


InertSearch™ for LC

Inertsil® Applications

Pionin

Data No. LA357-0000

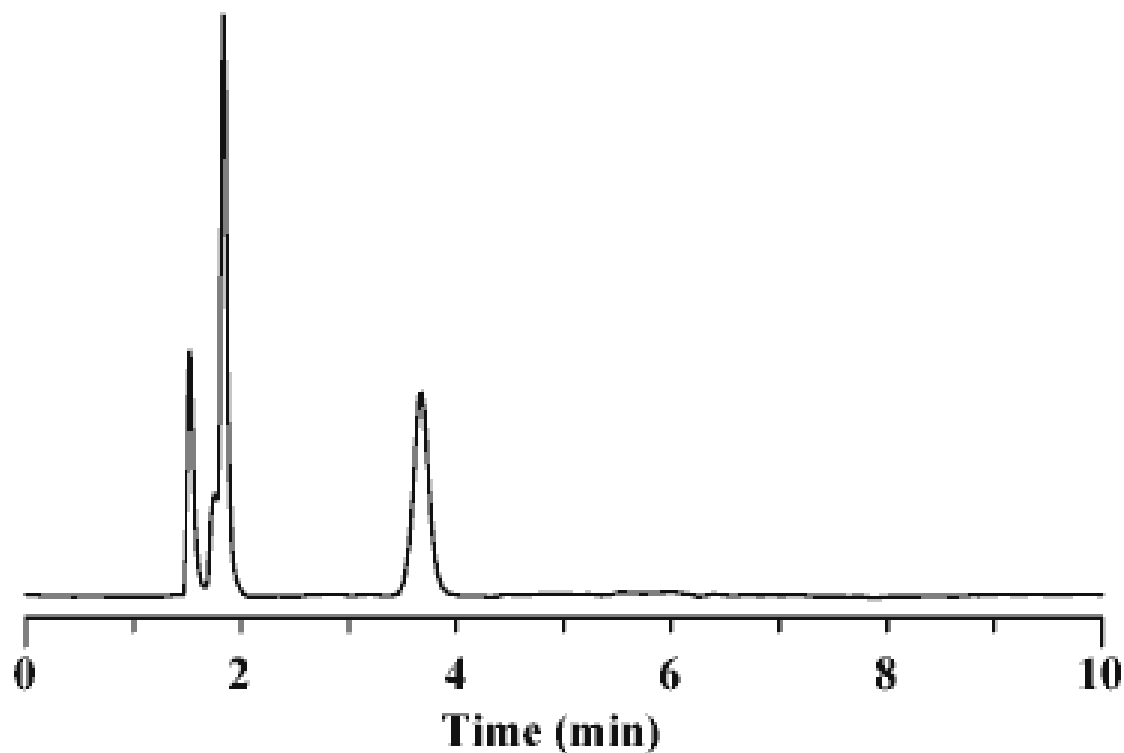

Column : Inertsil ODS-3V
150 × 4.6 mmI.D.

Solutes : 1. Pionin (0.04mg/mL)

Column Cat. No. :

Eluent : Methanol / Water (85 / 15) / acetic acid / triethylamine =
996 / 2 / 2

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254nm

Injection Vol. : 5µL

Sample : Pionin (0.04mg/mL)