


DNPH-Aldehydes

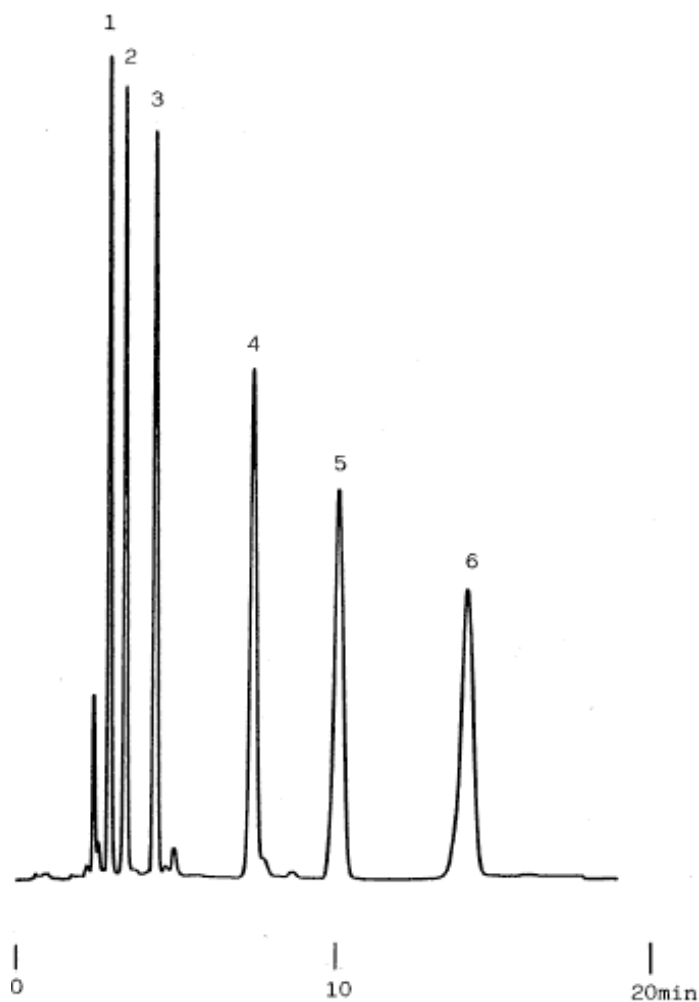

Column : Inertsil ODS-2
 5µm
 150 × 4.6 mm I.D.

Column Cat. No. : 5020-01101

Eluent : CH₃CN/H₂O = 70/30

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample : DNPH-Aldehyde

Solutes : 1. DNPH-Formaline (Formaldehyde)
 2. DNPH-Acetaldehyde
 3. DNPH-Propionaldehyde
 4. DNPH-n-Valeraldehyde
 5. DNPH-Capronaldehyde
 6. DNPH-n-Heptylaldehyde