


DNPH-Aldehydes

Data No. LA403-0000

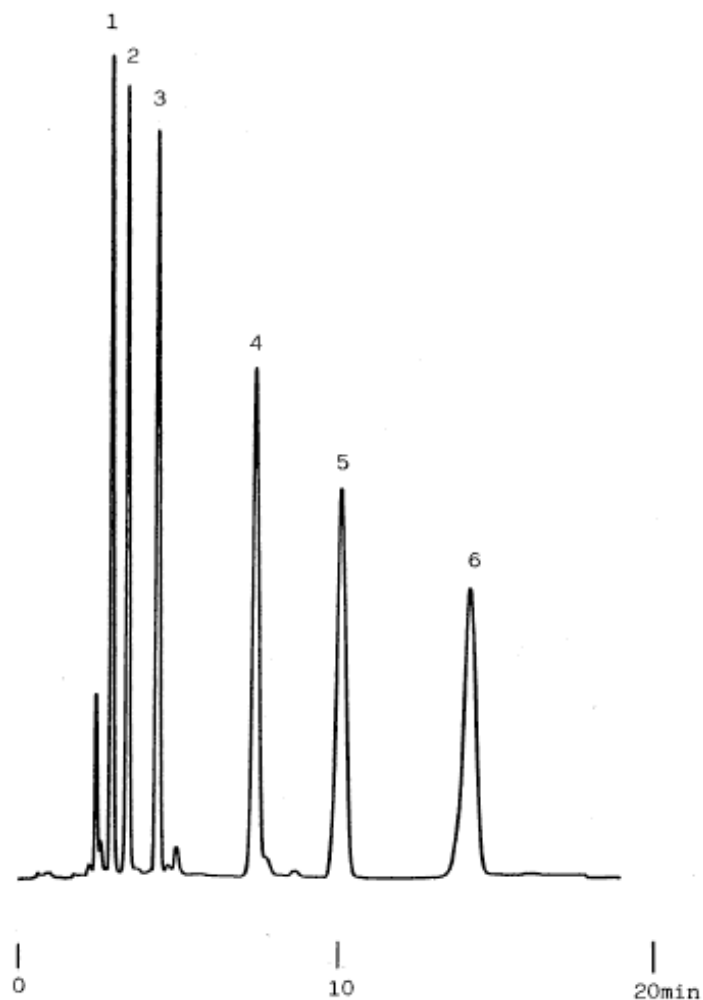

Column : Inertsil ODS-2
5 μ m
150 \times 4.6 mm I.D.

Column Cat. No. : 5020-01101

Eluent : CH₃CN/H₂O = 70/30

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample : DNPH-Aldehyde

Solutes : 1. DNPH-Formaline (Formaldehyde)
2. DNPH-Acetaldehyde
3. DNPH-Propionaldehyde
4. DNPH-n-Valeraldehyde
5. DNPH-Capronaldehyde
6. DNPH-n-Hepthylaldehyde