


Proteins

Data No. LA407-0000

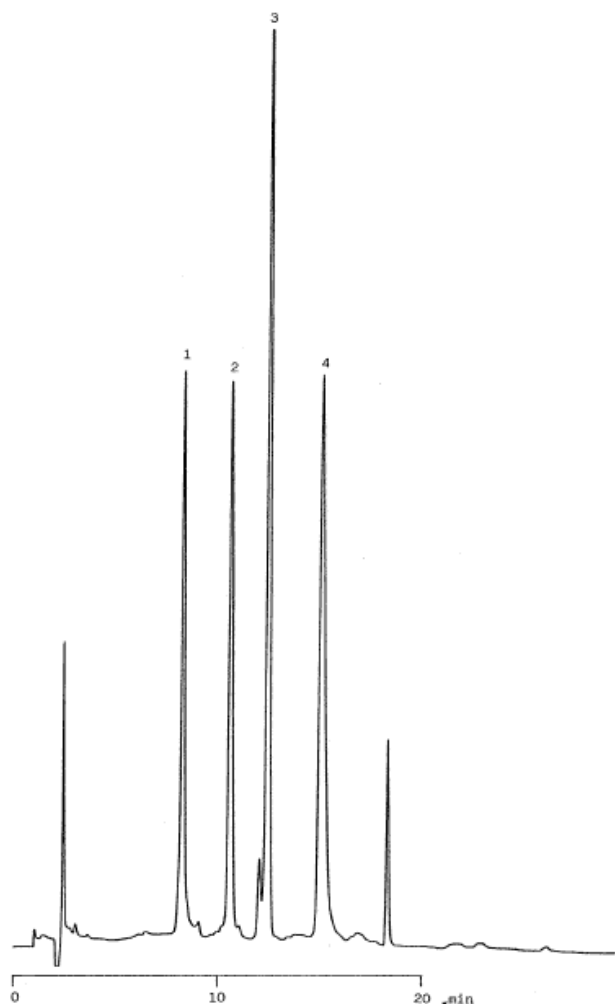

Column : Inertsil Ph
5 μ m
150 \times 4.6 mm I.D.

Column Cat. No. : 5020-01301

Eluent : A: CH₃CN
B: 0.1 % TFA
A/B = 20/80 - 20 min. - 60/40

Flow rate : 1.0 mL/min.

Col. Temp. : Ambient

Detection : UV 220 nm

Injection Vol. :

Sample : Protein

Solutes : 1. Ribonuclease A
2. Cytochrome C
3. Lysozyme
4. Myoglobin