


Tryptic Digest of BSA

Data No. LA410-0000

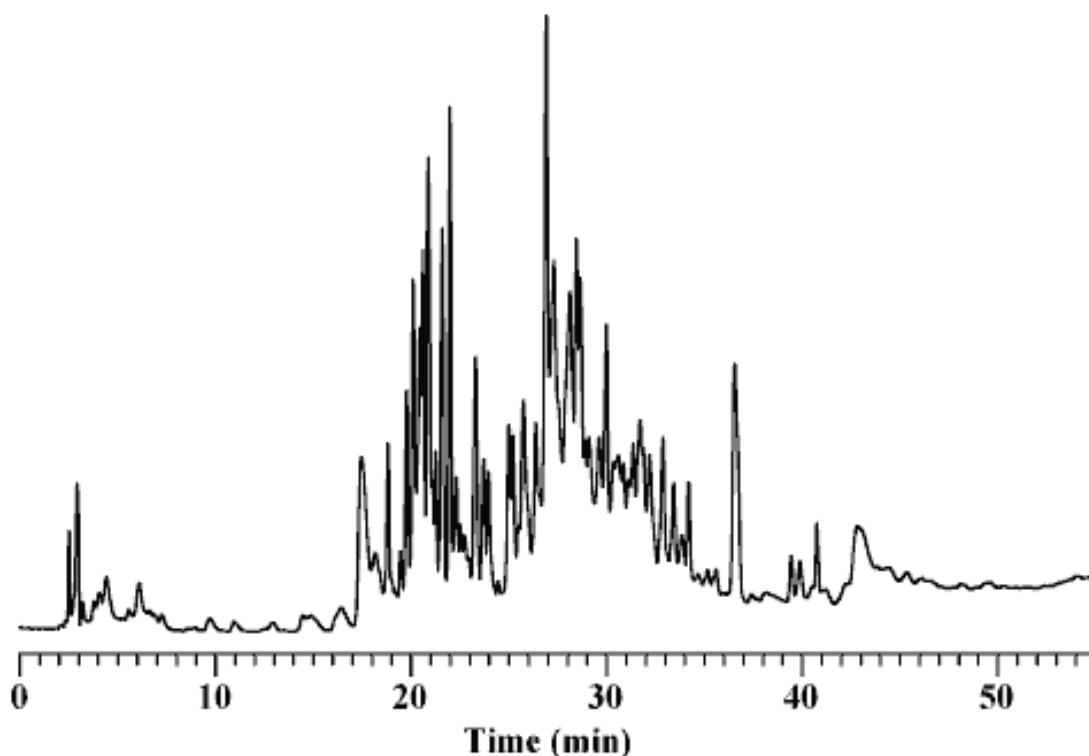

Column : Inertsil ODS-3
3 μ m
150 \times 0.3 mm I.D.

Solutes : 1. Tryptic Digest of BSA (5 mg/mL, 37°C, 24 h)

Column Cat. No. : 5020-01871

Eluent : A) 0.1 % TFA in CH₃CN/H₂O = 2/98
B) 0.1 % TFA in CH₃CN/H₂O = 90/10
A/B = 90/10 - 5 min. - 90/10 - 50 min. - 40/60

Flow rate : 5 μ L/min.

Col. Temp. : 30°C

Detection : UV 214 nm

Injection Vol. :

Sample : Tryptic Digest of BSA