


Anions

Data No. LA413-0000

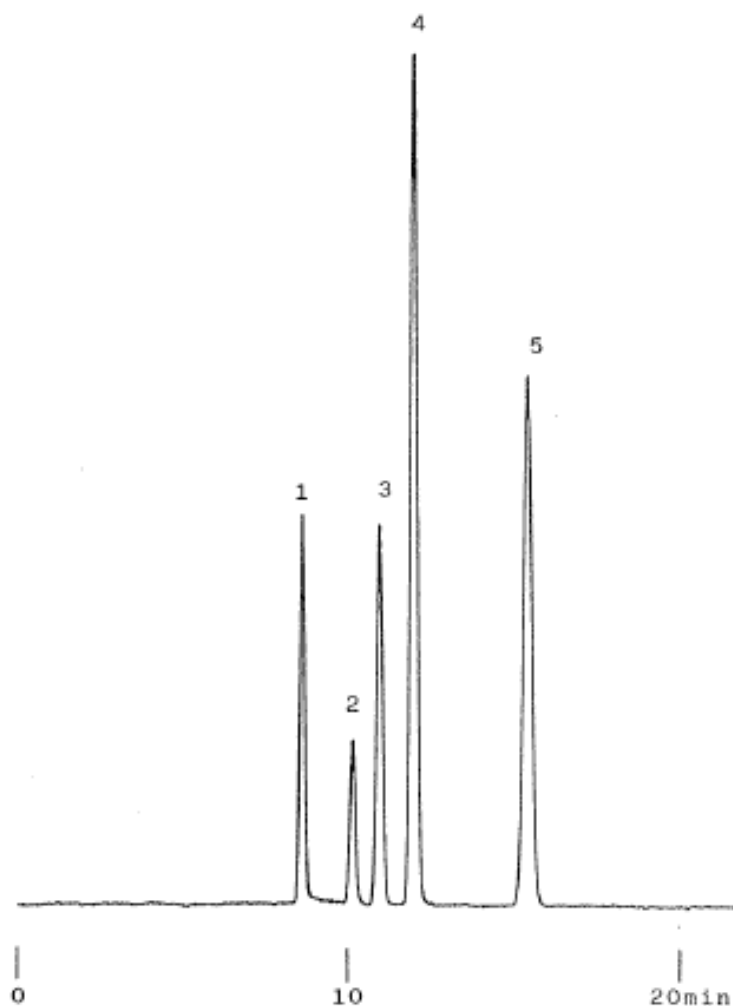

Column : Inertsil ODS-2
5 μ m
250 \times 4.6 mm I.D.

Column Cat. No. : 5020-01102

Eluent : CH₃CN/2 mM IPCA-03H (pH 7.0; HCl) = 20/80
IPCA-03H: Tetra-n-butylammonium Hydroxide

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 215 nm

Injection Vol. :

Sample : Inorganic anion

Solutes : 1. IO₃⁻
2. BrO₄⁻
3. Br⁻
4. NO₃⁻
5. I⁻