


Metal complexes of Etioporphyrin

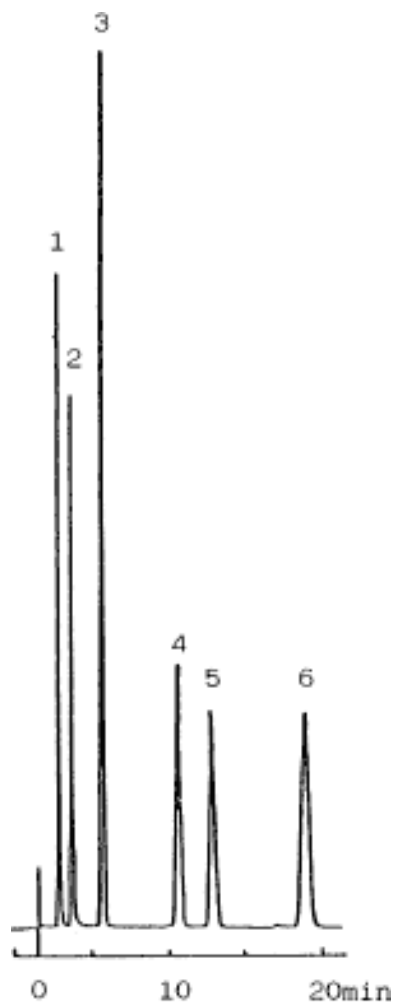

Column : Inertsil ODS-2
 5µm
 150 × 4.6 mm I.D.

Column Cat. No. : 5020-01101

Eluent : CH₃OH/Acetone = 80/20

Flow rate : 1.0 mL/min.

Col. Temp. : Ambient

Detection : VIS 400 nm

Injection Vol. :

Sample : Metal complex of etioporphyrin

Solutes : 1. Zn (etio)
 2. Vo (etio)
 3. H₂ (etio)
 4. Ni (etio)
 5. Pb (etio)
 6. Cu (etio)