


Pyrones

Data No. LA417-0000

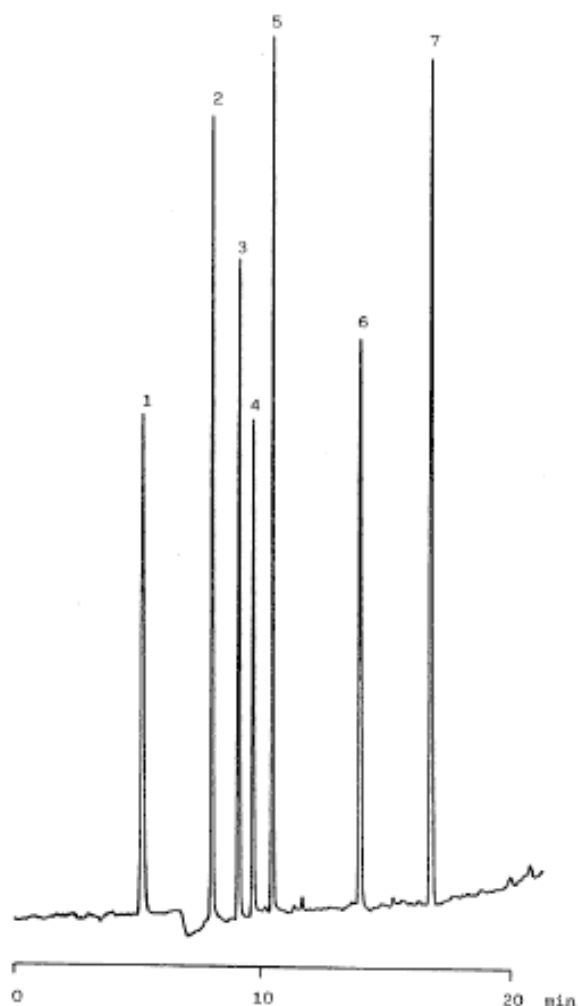

Column : Inertsil C4
5µm
150 × 4.6 mm I.D.

Column Cat. No. : 5020-01401

Eluent : A: CH₃CN
B: 0.05 % H₃PO₄
A/B = 0/100 - 20 min - 50/50

Flow rate : 1.0 mL/min.

Col. Temp. :

Detection : UV 270 nm

Injection Vol. :

Sample : Pyrone

Solutes : 1. 5-Hydroxy-2-hydroxymethyl-4-pyrone
2. 4-Pyrone-2-carboxylic acid
3. 2-Pyrone-5-carboxylic acid
4. 3-Hydroxy-2-methyl-4-pyrone
5. 4-Hydroxy-6-methyl-2-pyrone
6. 3-Hydroxy-2-methyl-4-pyrone propionate
7. 3-Hydroxy-2-methyl-4-pyrone isobutyrate