


Purification of C60 Buckminsterfullerene

Data No. LA420-0000

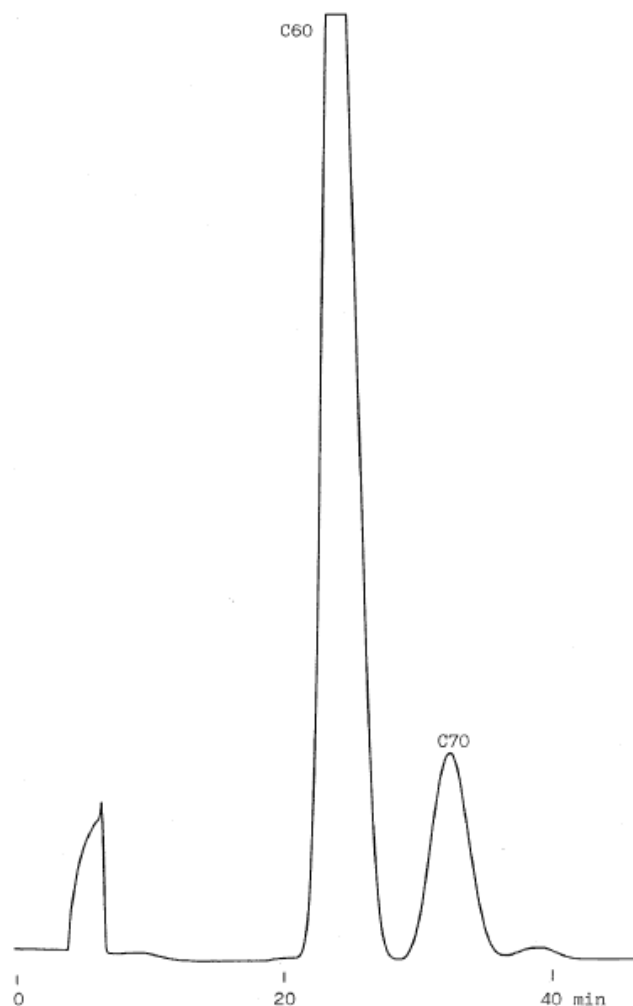

Column : ECONOPREP ODS
40 μ m
500 \times 100 mm I.D.

Solutes : 1. C60 (0.01 % in Eluent)
2. C70 (0.01 % in Eluent)

Column Cat. No. : 5020-36085

Eluent : Toluene/CH₃CN = 60/40

Flow rate : 280 mL/min.

Col. Temp. :

Detection : UV 350 nm

Injection Vol. : 1000 μ L

Sample : Crude C60, C70