


## Anions

Data No. LA421-0000

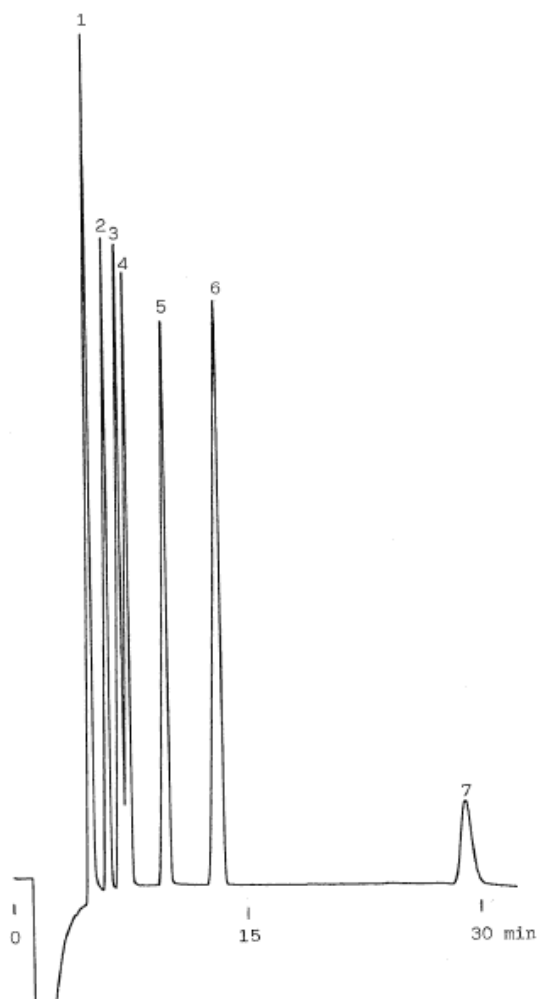

Column : Inertsil C4  
5 $\mu$ m  
150  $\times$  4.6 mm I.D.

Column Cat. No. : 5020-01401

Eluent : 1 mM IPCA-03H + 0.5 mM Phthalic acid

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample : Inorganic anion

Solutes : 1. F-  
2. Cl-  
3. Br-  
4. NO<sub>2</sub>-  
5. NO<sub>3</sub>-  
6. ClO<sub>3</sub>-  
7. SO<sub>4</sub><sup>2-</sup>