


Anions in Rain water

Data No. LA423-0000

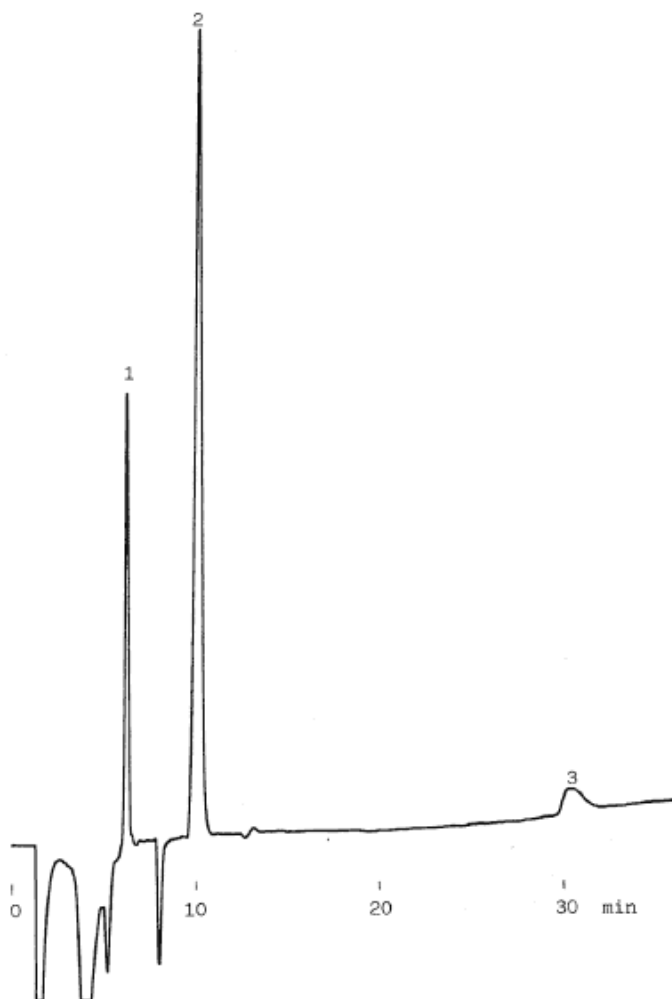

Column : Inertsil C4
5 μ m
150 \times 4.6 mm I.D.

Solutes : 1. Cl⁻
2. NO₃⁻
3. SO₄²⁻

Column Cat. No. : 5020-01401

Eluent : 1 mM IPCA-03H + 0.5 mM Phthalic acid

IPCA-03H: Tetra-n-butylammonium Hydroxide

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample : Inorganic anion