


Divinylbenzene

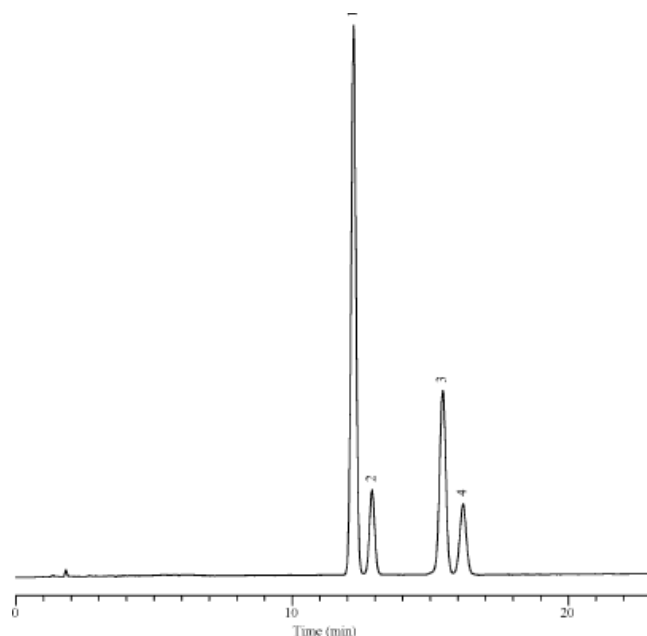

No.	R.Time	Area	Area %
1	12.230	423533	56.7
2	12.896	69140	9.3
3	15.442	181250	24.3
4	16.191	73337	9.8

Column : Inertsil ODS-P
5µm
250 × 4.6 mm I.D.

Solutes : Divinylbenzene Mixture

Column Cat. No. : 5020-02002

Eluent : CH₃CN/H₂O = 65/35

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample : Commercially available divinylbenzene