


## Analysis of Benzene Derivatives

Data No. LA431-0000

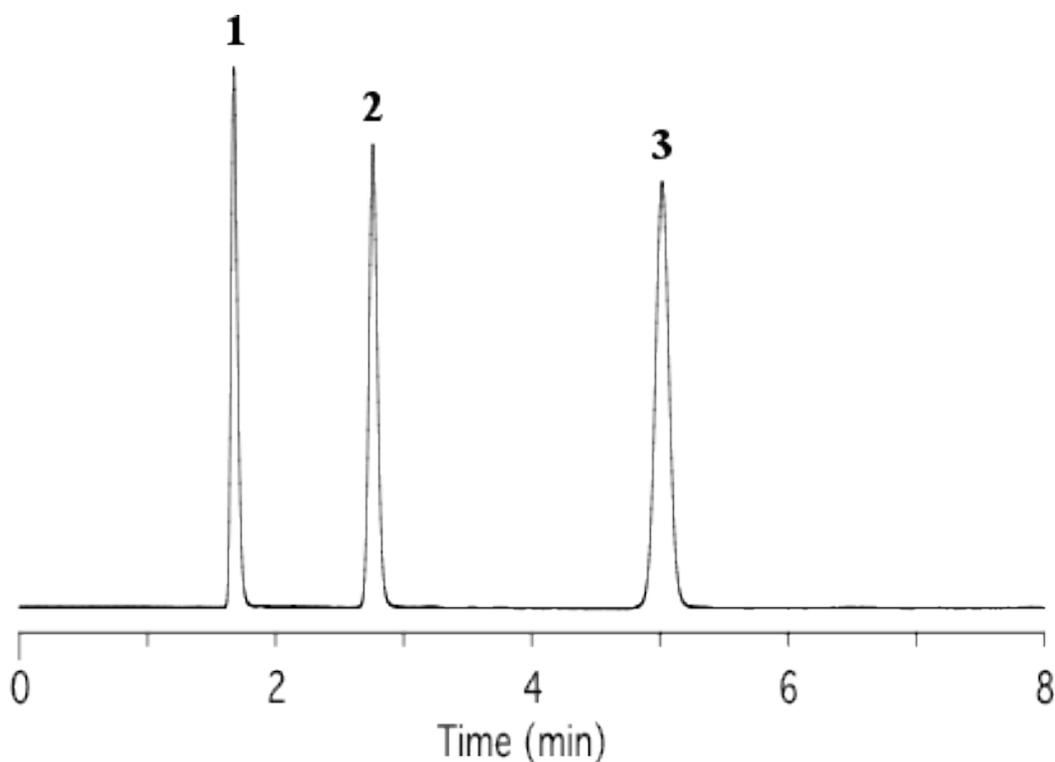

Column : Inertsil Diol  
5  $\mu$ m  
150 x 4.6 mm I.D.

Solutes : 1. Benzene  
2. Nitrobenzene  
3. o-Nitroanisole

Column Cat. No. : 5020-05645

Eluent : A) n-Hexane  
B) Ethanol  
A/B = 95/5, w/w

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample : Analysis of Benzene Derivatives