


Analysis of Benzene Derivatives

Data No. LA434-0000

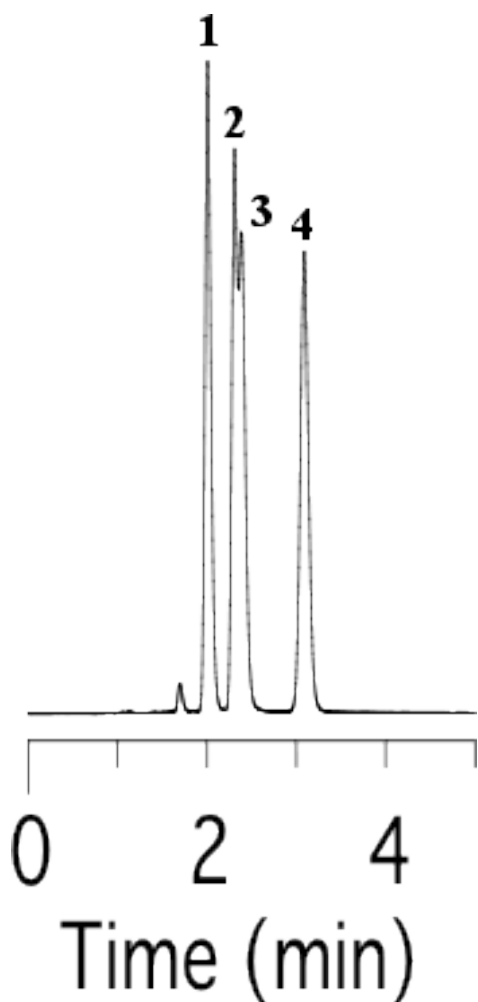

Column : Inertsil Diol
5 μ m
150 x 4.6 mm I.D.

Column Cat. No. : 5020-05645

Eluent : A) CH₃CN
B) H₂O
A/B = 30/70, w/w

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes : 1. Acetophenone
2. Benzene
3. Toluene
4. Naphthalene