


## Analysis of Aromatic Compounds

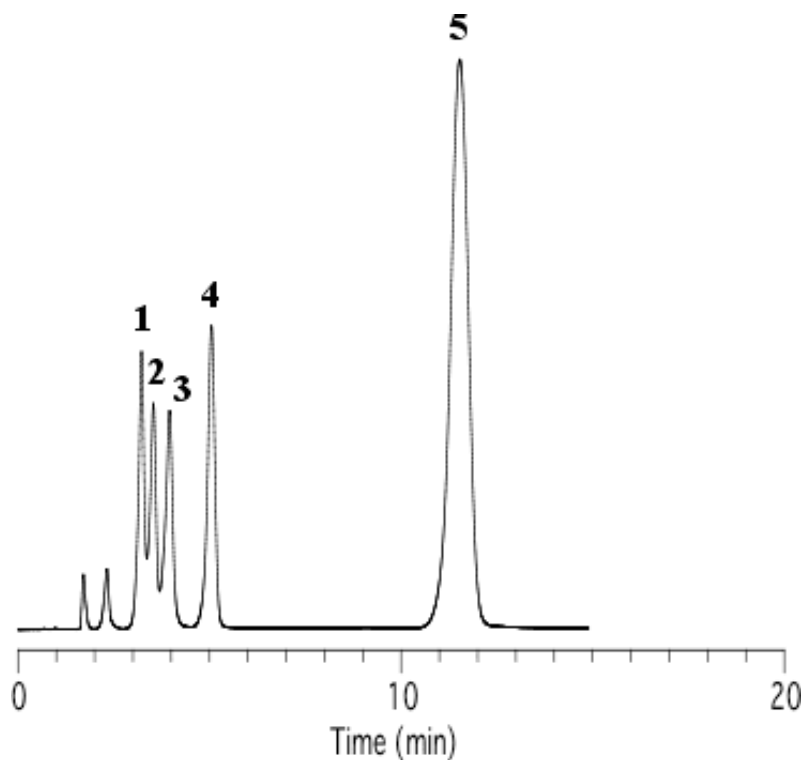

Column	: Inertsil Diol 5 $\mu$ m 150 x 4.6 mm I.D.	Solutes	: 1. Ethylbenzene 2. n-Propylbenzene 3. n-Butylbenzene 4. Naphthalene 5. Anthracene
Column Cat. No.	: 5020-05645		
Eluent	: A) CH <sub>3</sub> CN B) H <sub>2</sub> O A/B = 40/60, v/v		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	:		
Sample	:		