


Analysis of Samples for Column Selectivity

Data No. LA441-0000

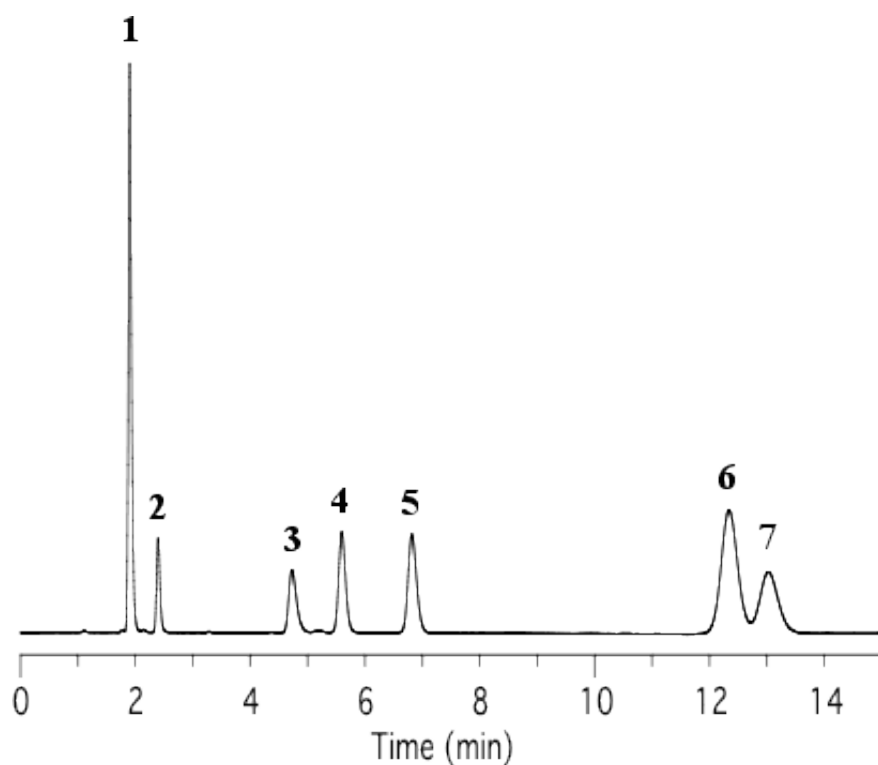

Column	: Inertsil Ph-3 5 µm 150 x 4.6 mm I.D.	Solutes	: 1. Uracil 2. Phenol 3. Caffeine 4. n-Butylbenzene 5. n-Amylbenzene 6. Triphenylene 7. o-Terphenyl
Column Cat. No.	: 5020-01302		
Eluent	: A) CH ₃ OH B) H ₂ O A/B = 60/40, v/v		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	:		
Sample	:		