


Analysis of Aromatic Compounds

Data No. LA445-0000

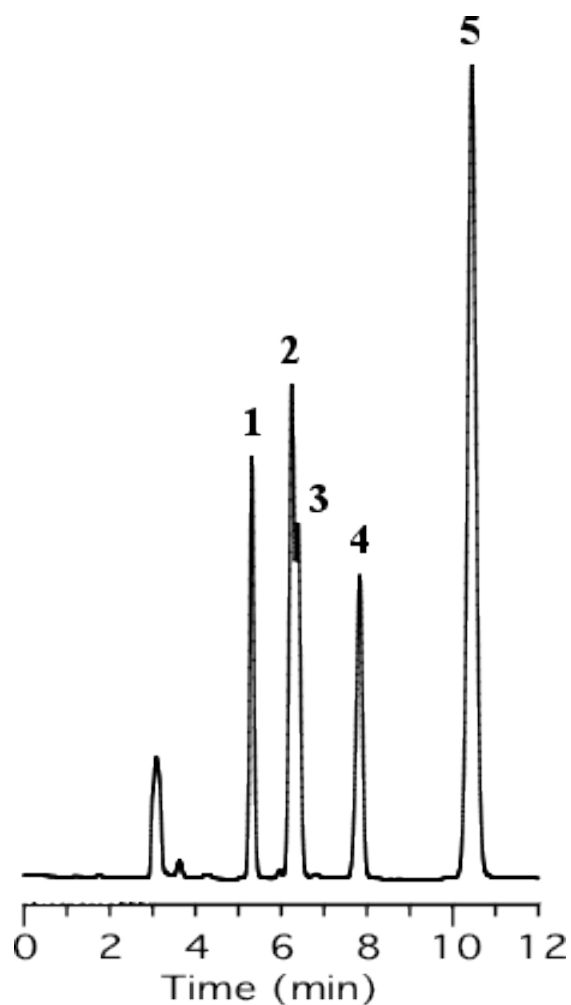

Column : Inertsil CN-3
5 μ m
150 x 4.6 mm I.D.

Column Cat. No. : 5020-01940

Eluent : A) CH₃CN
B) H₂O
A/B = 40/60, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes : 1. Ethylbenzene
2. Naphthalene
3. n-Propylbenzene
4. n-Butylbenzene
5. Anthracene