


Analysis of Benzene Derivatives

Data No. LA460-0000

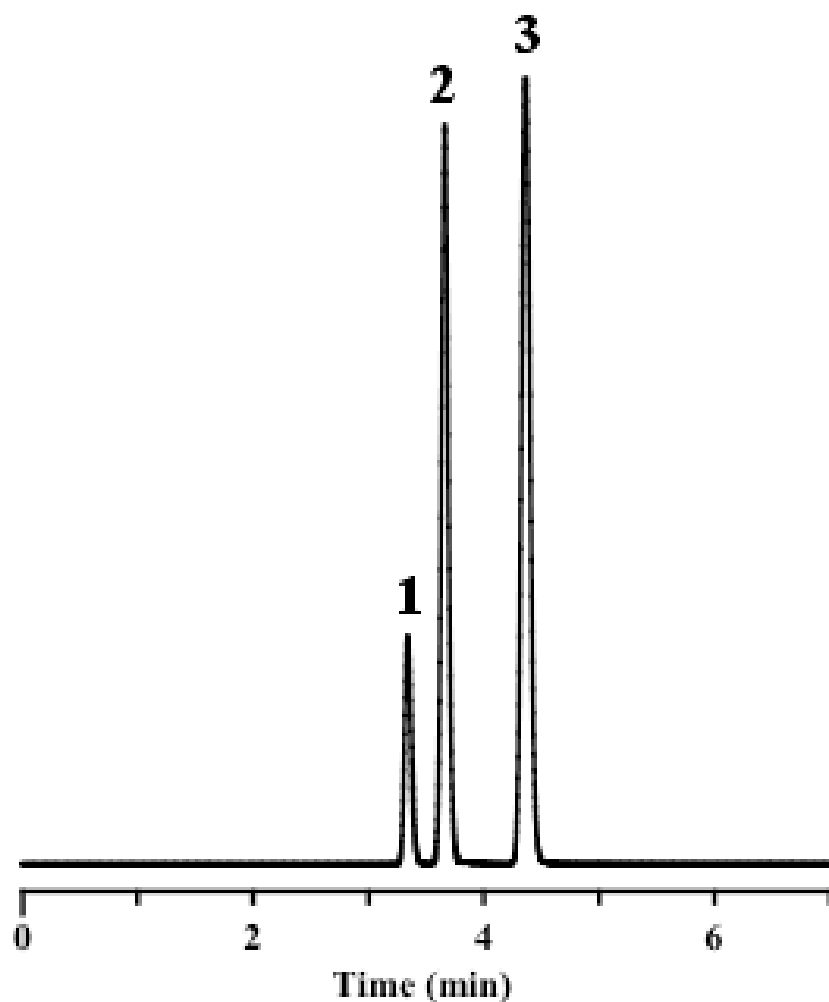

Column : Inertsil WP300 SIL
5 μm
250 x 4.6 mm I.D.

Solutes : 1. Benzene
2. Nitrobenzene
3. o-Nitroanisole

Column Cat. No. : 5020-06046

Eluent : A) n-Hexane
B) Ethanol
A/B = 98/2, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :