


Analysis of Divinylbenzene (Mixed Isomers)

Data No. LA475-0000

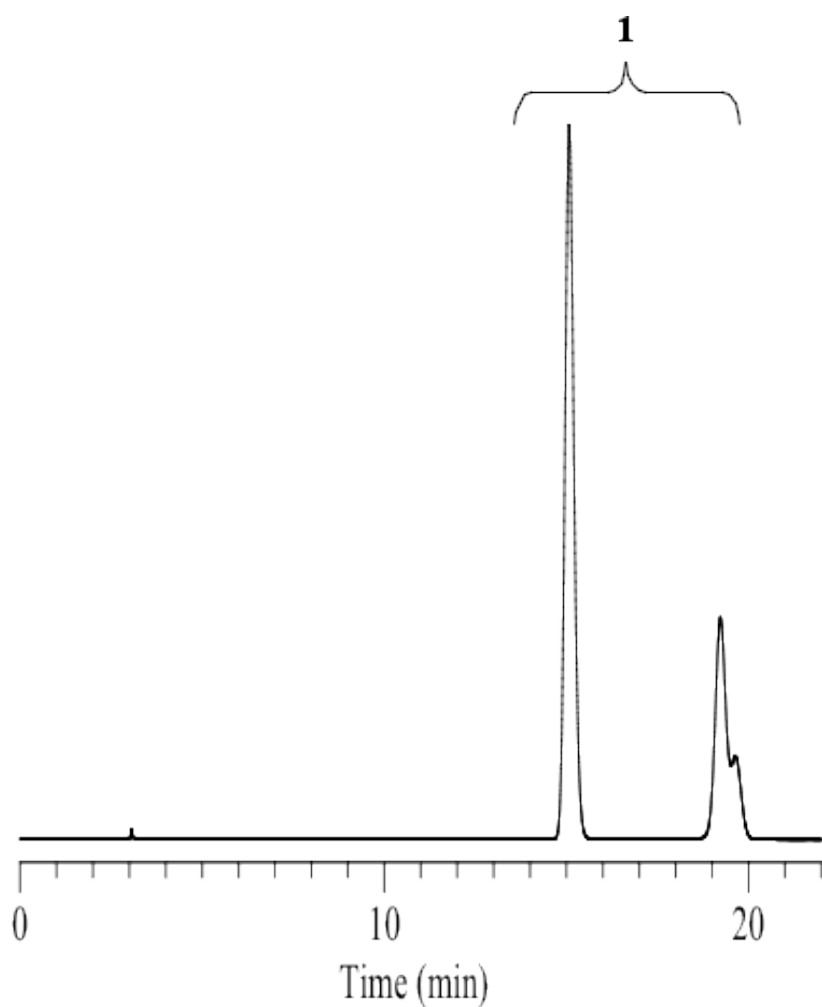

Column : Inertsil ODS-3
5 μ m
250 x 4.6 mm I.D.

Column Cat. No. : 5020-01732

Eluent : A) CH₃CN
B) H₂O
A/B = 65/35, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes : 1. Divinylbenzene (Mixed Isomers)