


Analysis of Xylene Isomers

Data No. LA476-0000

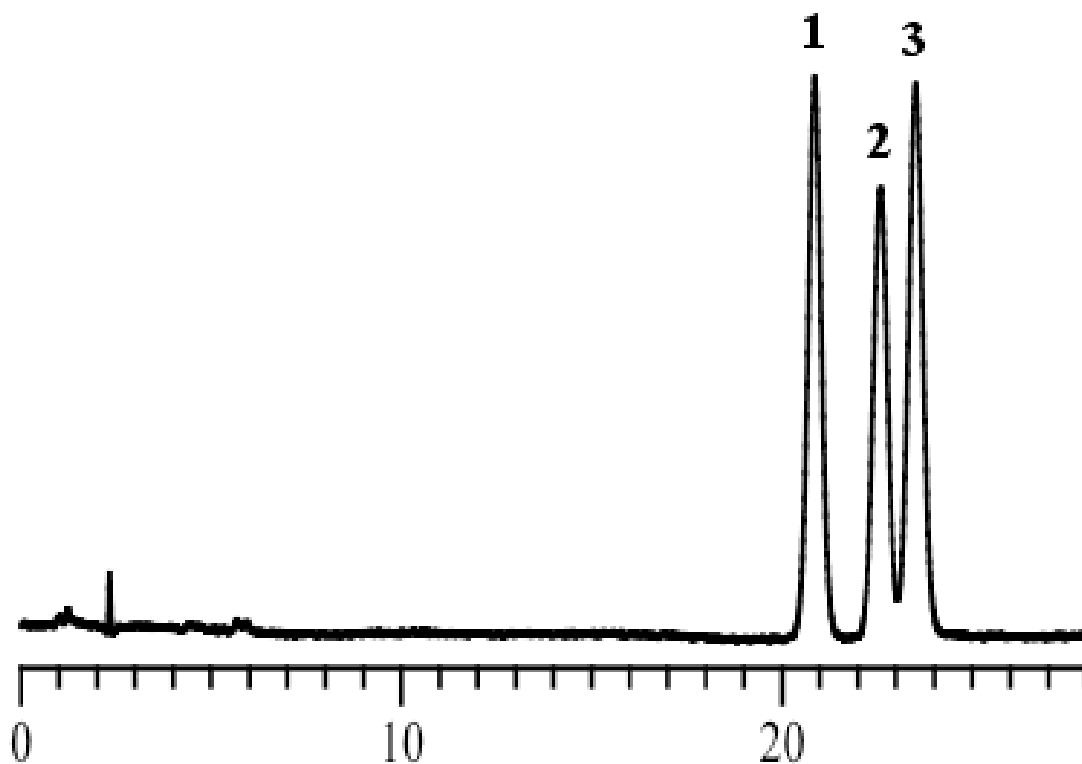

Column	: Inertsil ODS-P 5 µm 250 x 4.6 mm I.D.	Solutes	: 1. o-Xylene 2. m-Xylene 3. p-Xylene
Column Cat. No.	: 5020-02002		
Eluent	: A) Ethanol B) β-Cyclodextrin/H ₂ O = 5/995, v/v A/B = 50/50, v/v		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	:		
Sample	:		