


Analysis of Pyridine

Data No. LA486-0000

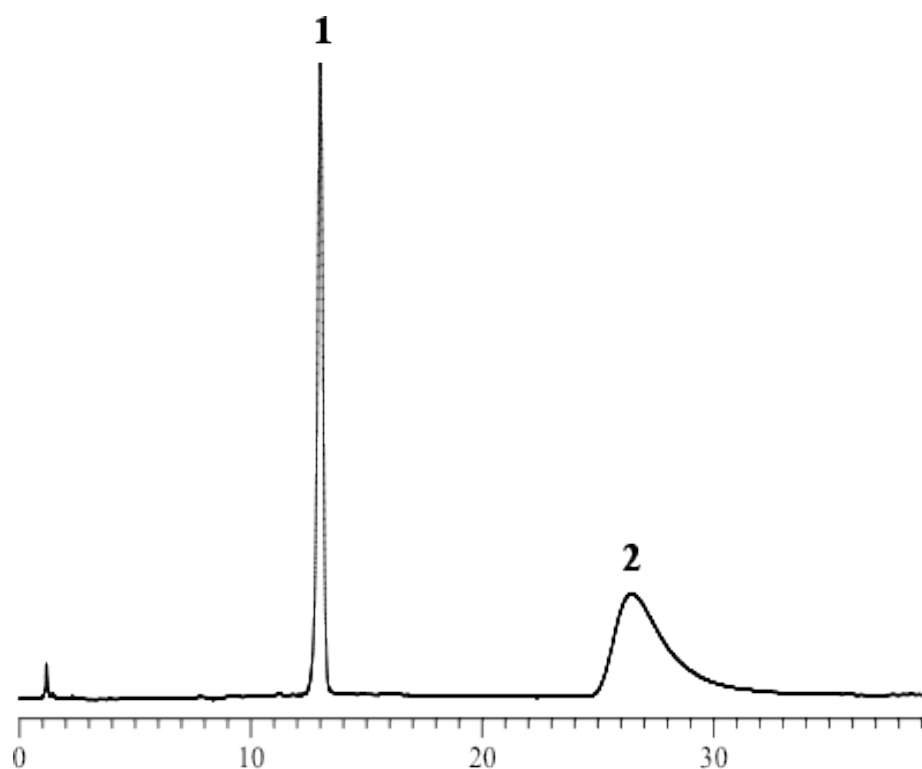

Column	: Inertsil ODS-P 5 µm 250 x 4.6 mm I.D.	Solutes	: 1. Phenol 2. Pyridine
Column Cat. No.	: 5020-02002		
Eluent	: A) CH ₃ OH B) H ₂ O A/B = 30/70, v/v		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	:		
Sample	:		