


# Analysis of Aromatic Compounds

Data No. LA487-0000

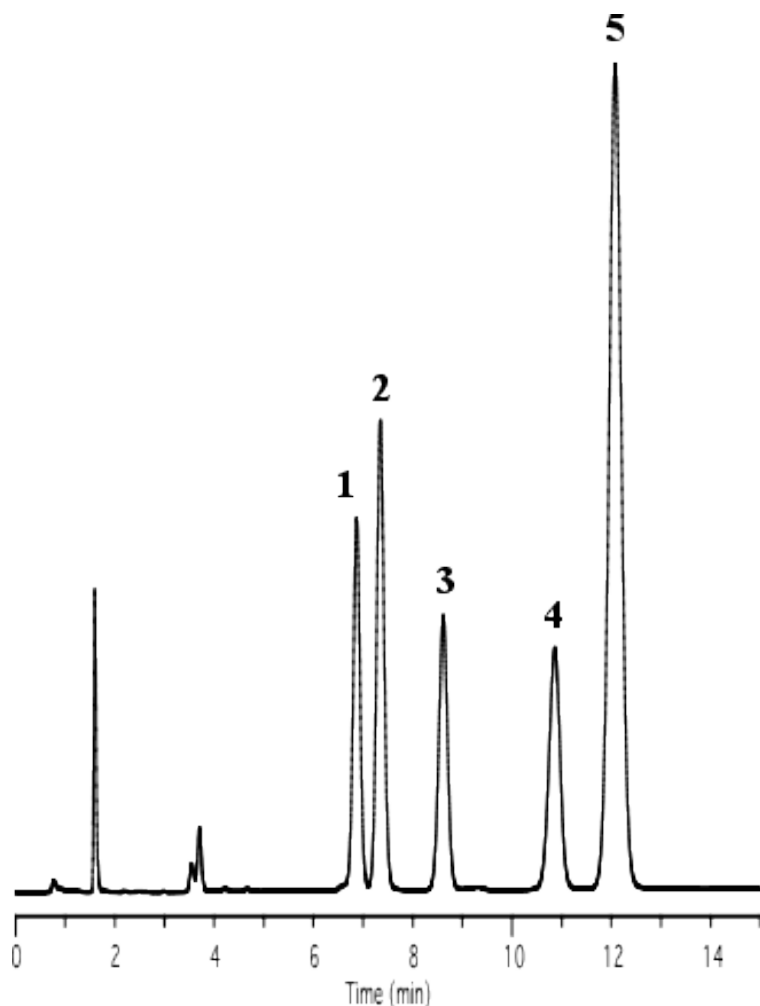

Column : Inertsil Ph-3  
5  $\mu$ m  
150 x 4.6 mm I.D.

Column Cat. No. : 5020-01302

Eluent : A) CH<sub>3</sub>CN  
B) H<sub>2</sub>O  
A/B = 50/50, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes : 1. Ethylbenzene  
2. Naphthalene  
3. n-Propylbenzene  
4. n-Butylbenzene  
5. Anthracene