


Analysis of Aromatic Compounds

Data No. LA488-0000

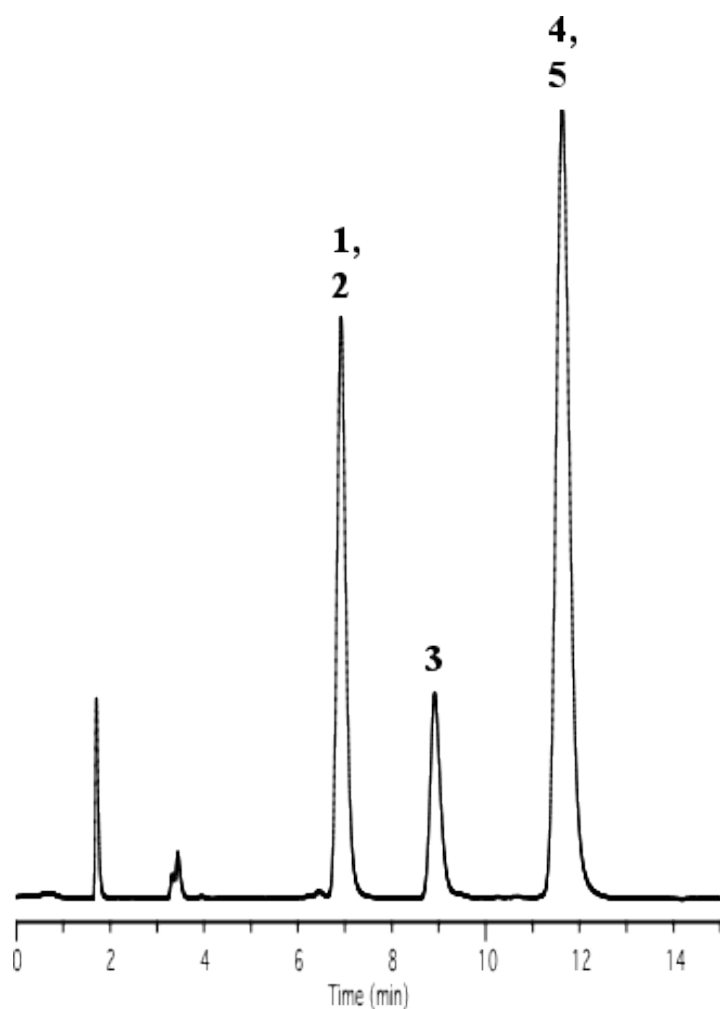

Column : Inertsil Ph
5 μ m
150 x 4.6 mm I.D.

Column Cat. No. : 5020-01301

Eluent : A) CH₃CN
B) H₂O
A/B = 50/50, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes : 1. Ethylbenzene
2. Naphthalene
3. n-Propylbenzene
4. n-Butylbenzene
5. Anthracene