


InertSearch for LC Data No. LA488-0000

Analysis of Aromatic Compounds

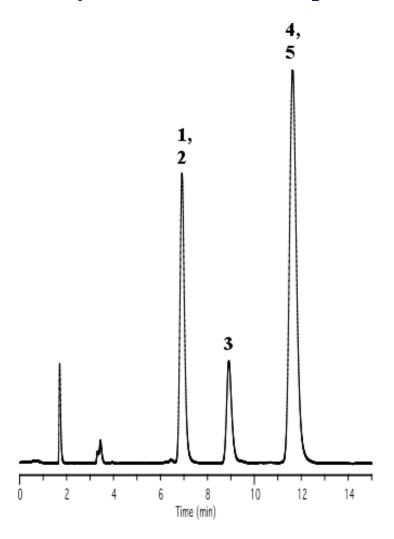

Column : Inertsil Ph

5 µm

150 x 4.6 mm I.D.

Column Cat. No.: 5020-01301

Eluent : A) CH₃CN

B) H₂O

A/B = 50/50, v/v

Flow rate :1.0 mL/min.

Col. Temp. :40°C

Detection : UV 254 nm

Injection Vol. :
Sample :

Solutes

_

2. Naphthalene

: 1. Ethylbenzene

3. n-Propylbenzene

4. n-Butylbenzene

5. Anthracene