


Analysis of α -Lipoic acid

Data No. LA512-0696

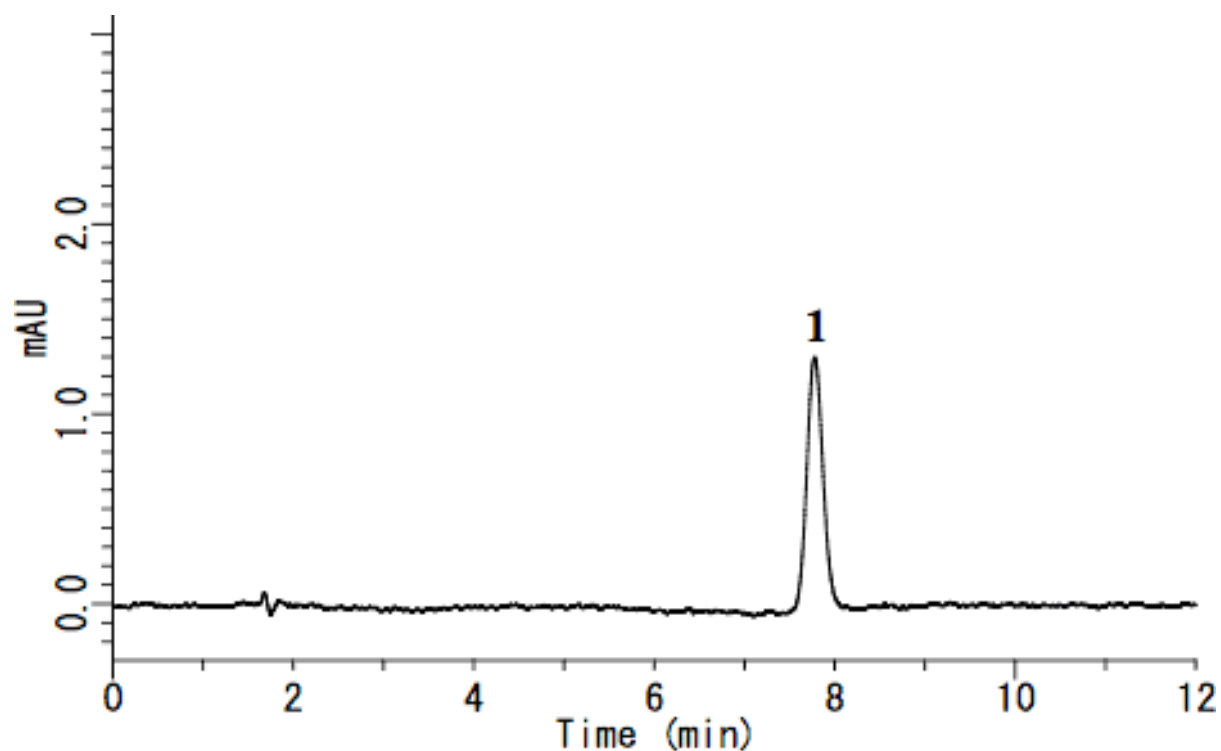

System : GL-7400 system

Solutes : 1. α -Lipoic acid (1.0 mg/L)

Column : Inertsil ODS-3 (4 μ m, 150 x 4.6 mm I.D.)

Column Cat. No. :

Eluent : A) CH₃OH

B) 20 mM NaH₂PO₄ (pH 2.5; H₃PO₄)

A/B = 60/40 ,v/v

Flow rate : 1.0 mL/min

Col. Temp. : 40 °C

Detection : UV 333 nm (GL-7450 UV Detector)

Injection Vol. : 20 μ L

Sample :