


## Analysis of Bufadienolides

Data No. LA544-0698

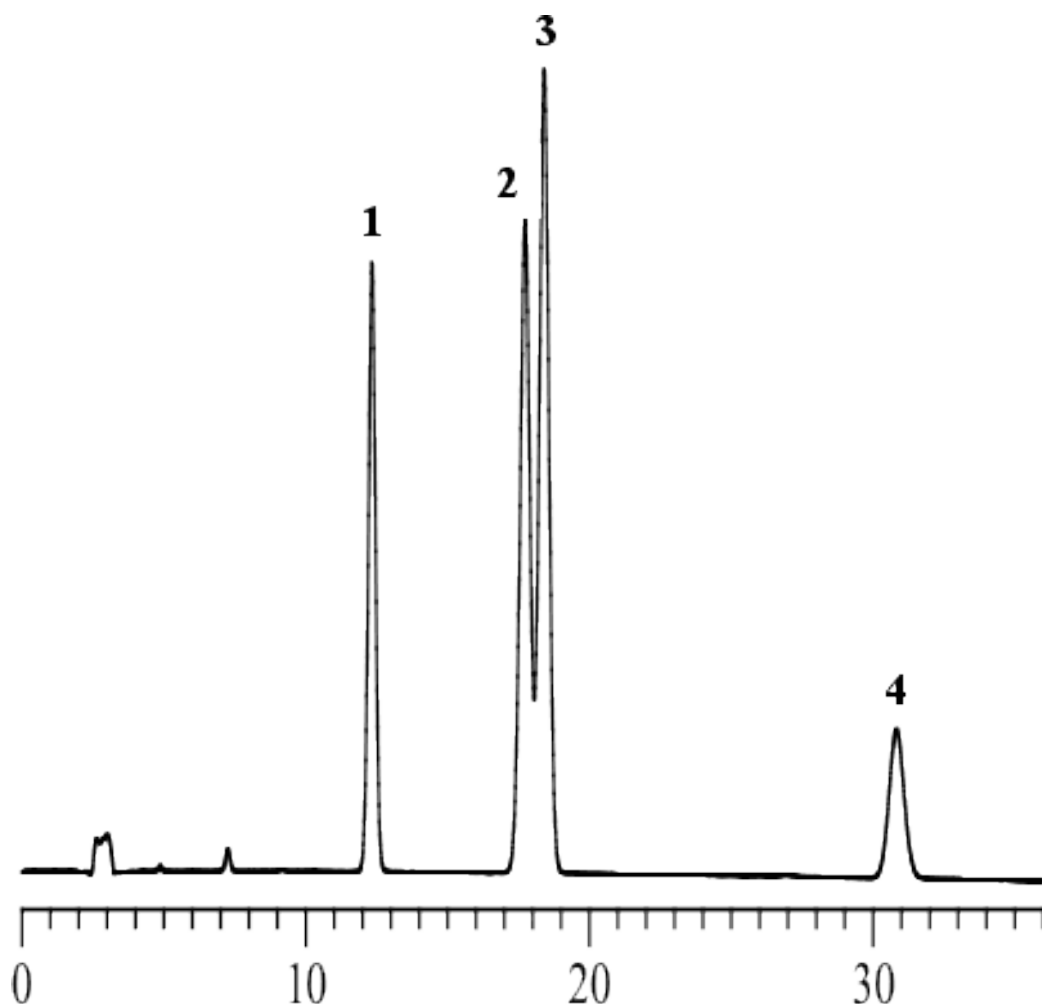

Column : Inertsil ODS-3  
5 µm  
250 x 4.6 mm I.D.

Column Cat. No. :

Eluent : A) CH<sub>3</sub>CN  
B) 0.1 % H<sub>3</sub>PO<sub>4</sub>  
A/B = 45/55, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 300 nm

Injection Vol. :

Sample :

Solutes : 1. Bufalin  
2. Cinobufagin  
3. Resibufogenin  
4. Indomethacin