


InertSearch™ for LC

Inertsil® Applications

Analysis of Protein (Inertsil WP300 Diol)

Data No. LA672-0699

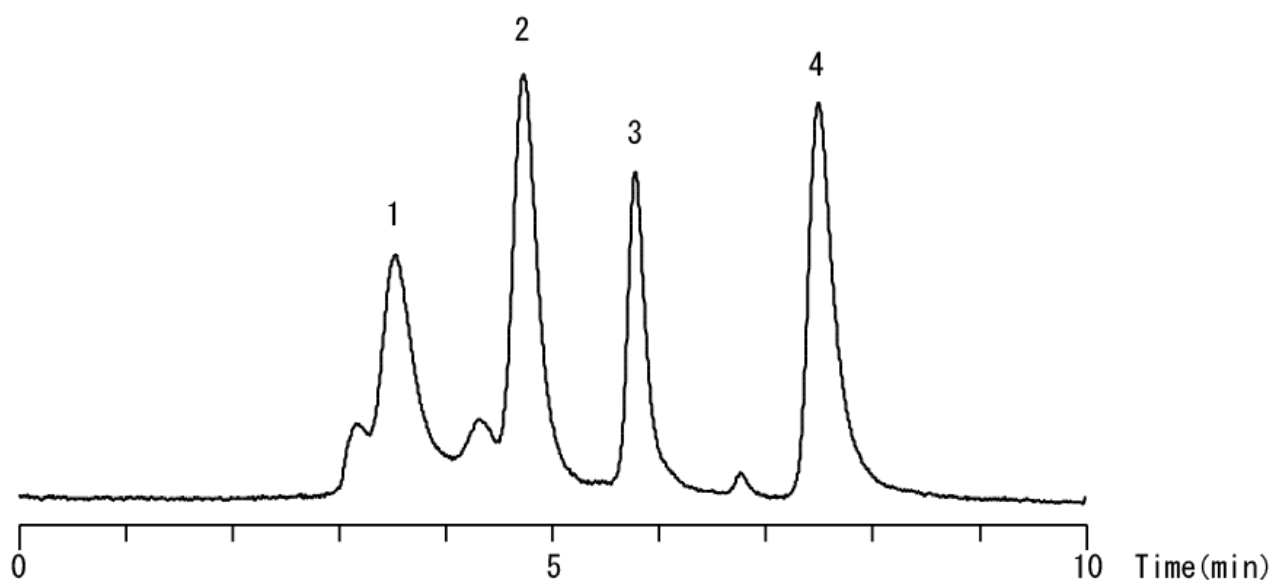

Conditions

Column : Inertsil WP300 Diol (5 µm, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-05746
Eluent : 0.1 M Na₂HPO₄ (pH 6.9, NaH₂PO₄)
Flow Rate : 0.5 mL/min
Col. Temp. : 30 °C
Detection : UV 220 nm
Injection Vol. : 20 µL

Analyte

1. Thyrogloblin (0.25 mg/mL)
2. BSA (0.25 mg/mL)
3. Insulin Chain A (0.25 mg/mL)
4. Oxytocin (0.21 mg/mL)