


# InertSearch™ for LC

Inertsil® Applications

## Analysis of Vitamin K (Inertsil ODS-2)

Data No. LA765-0696

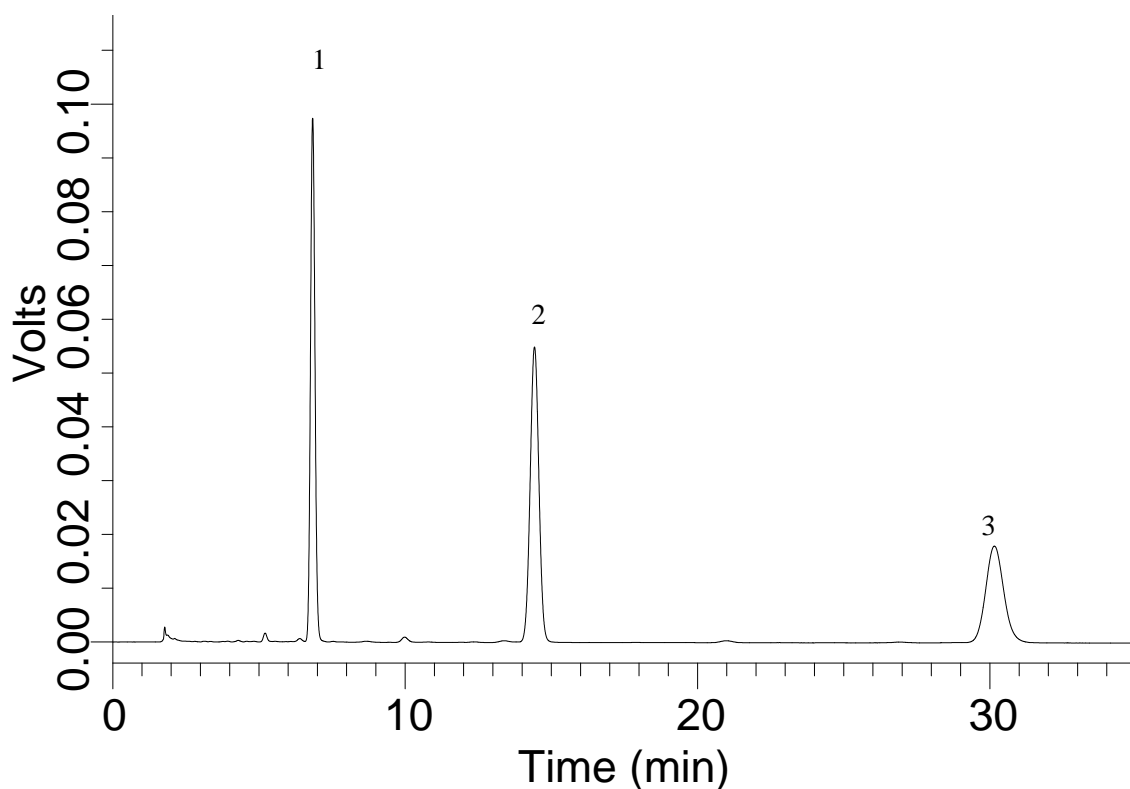

### Conditions

**System** : LC800 HPLC system  
**Column** : Inertsil ODS-2 (5  $\mu$  m, 150 x 4.6 mm I.D.)  
**Column Cat. No.** : 5020-01101  
**Eluent** : CH<sub>3</sub>CN  
**Flow Rate** : 1.0 mL/min  
**Col. Temp.** : 40 °C  
**Detection** : UV 270 nm (LC800 UV Detector)

### Analyte

1. Vitamin K2(MK-4)  
2. Vitamin K1  
3. Vitamin K2(MK-7)