


InertSearch™ for LC

Inertsil® Applications

Analysis of Clioquinol (Inertsil ODS-4)

Data No. LA780-0521

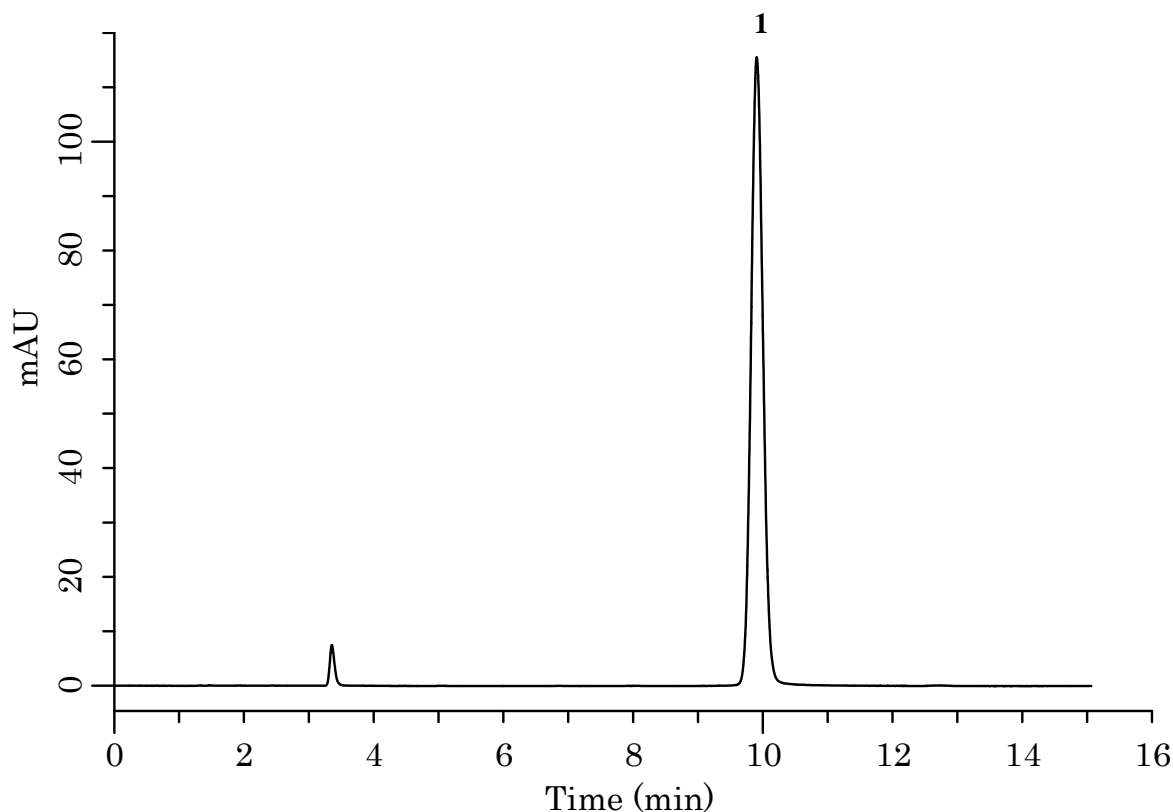

Conditions

System : GL-7400 HPLC system
Column : Inertsil ODS-4 (5 µm, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-03945
Eluent : A) CH₃CN
 : B) 0.1 % H₃PO₄
 : A/B = 60/40, v/v
Flow Rate : 1 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm (GL-7450 UV Detector)
Injection Vol : 2µL

Analyte

1. Clioquinol (100 µg/mL)