


InertSearch™ for LC

Inertsil® Applications

Analysis of Malachite green and Leucomalachite green

Data No. LA942-0696

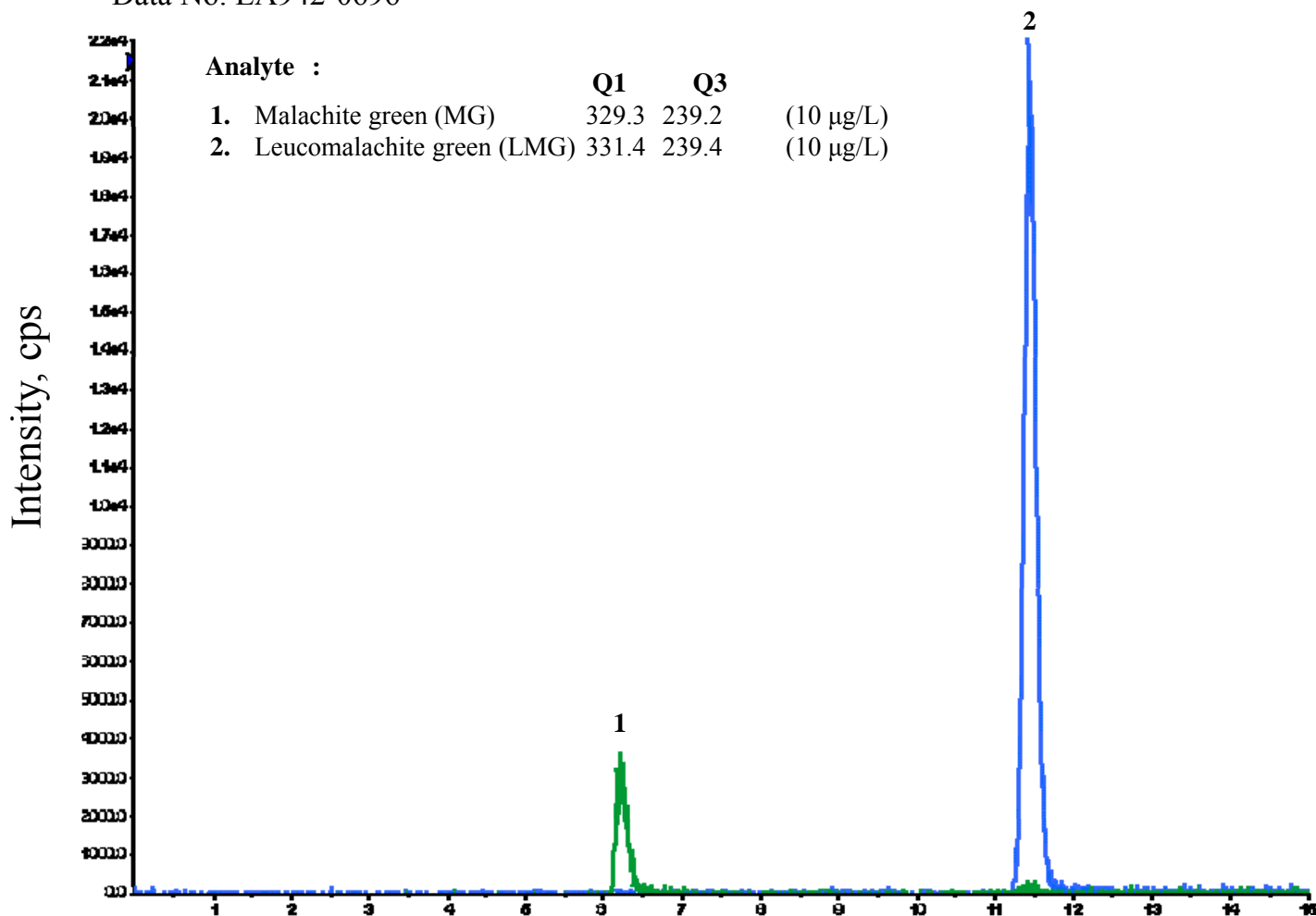

Conditions

System : Symbiosys Pharma system (LC mode)
API2000®

Column : Inertsil ODS-3 (3µm, 150 x 2.1 mm I.D.)

Column Cat. No. : 5020-04415

Eluent : A) CH₃CN
B) 0.1% HCOOH
A/B = 20/80 -10 min-95/5-10 min -95/5, v/v

Flow rate : 0.2 mL/min

Col. Temp. : 40 °C

Detection : LC/MS/MS (API2000®: ESI, Positive, MRM)

Injection Vol. : 10 µL

Sample : Standard